



HyperpureX®


FE

Eminent series

Medium laboratory water purification system



Inside Equals Outside Excellent Performance

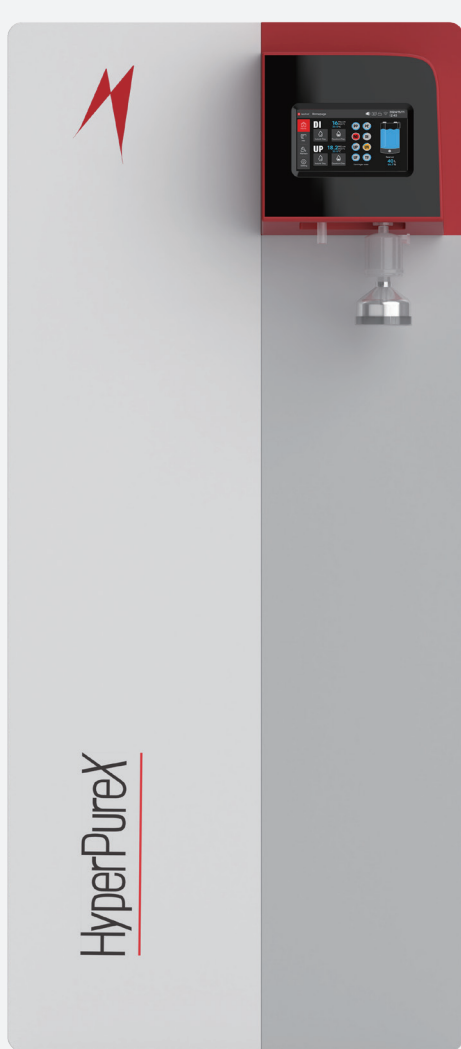

HyperpureX[®]

FE Eminent series (FEU/FED/FEUS/FEDS FEUE/FEDE/FERS)

Medium laboratory water purification system

Breakthrough design to highlight the aesthetics of science and technology.

FE series, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, rigorous double RO system^[1], advanced EDI module^[2] and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank, can meet your critical and professional application requirements of lab pure water.

System output ranges from 50 to 250 liters/h. It can simultaneously produce ultrapure water (18.2MΩ.cm), high pure water (>17.5MΩ.cm) or pure water^[3]. The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ASTM D1193-06, GB/T 11446.1-2013, GB/T 33087-2016, GB/T 6682-2008, CP, EP, USP, JP, CAP, CLSI, etc.

^[1] The double RO system is only used for FEUS/FEDS/FEUE/FEDE/FERS series products.

^[2] EDI module is only used for FEUE/FEDE series products.

^[3] FED/FERS series products can produce single RO water (ion rejection rate $\geq 98\%$). EDS/ERS series products can produce double RO water ($<5\mu\text{s/cm}$).

Application Area:

- HPLC、UPLC、LC-MS
- ICP-MS、ICP-AES、AAS、GC-MS
- MALDI-TOF-MS、IC、TOC analysis
- Electrochemical, spectrophotometric determination
- Preparation of microbial media and reagents
- Cell culture, PCR, IVF
- Protein purification, electrophoresis, biochemistry
- Proteomics, genomics, immunoassay
- Feed water of laboratory instruments, such as: autoclave, bottle washing machine, environmental test chamber, water bath, etc.

Innovative Control System Bringing Efficiency & Intelligence

Intelligent human-computer interactive control system

- Intelligent human-computer interactive control system
- 5-inch IPS colorful capacitive touch screen, resolution: 800×480, free switching in dual language (Chinese-English) display, to achieve touch operation experience the same as mobile phone.
- All new style intuitive UI interface design, convenient to fully understand the system operating status and parameters. A glance for all important information.

Traceable comprehensive data management

- Store operating data records up to 3 years^[1], including water dispensing, alarm and cartridge replacement, achieve data storage and paperless data management of the whole product life cycle by the cloud platform, meeting data tracking needs.
- With function of data exporting from USB port and data downloading from cloud platform.
- Data report documents in EXCEL format, meet traceability provisions of data integrity.
- Function of water dispensing report, including each water quality, water volume and user information, meet the regulatory requirements, easier to certificate.

^[1] Storage duration varies depending on dispensing and alarm frequency.

Innovative Internet of Things (IOT) and cloud platform technology

- Access the internet by 4G or WIFI, achieve remote data acquisition, monitoring and management.
- Log into the cloud platform from PC, WAP or WeChat to get the device information.
- Timely alarm information to achieve fast customer service response.
- Health analysis based on big data makes fault judgment more accurate.
- Work order system can provide customers with efficient service guarantee.
- It can be connected to LIMS or BMS to realize equipment informatization and make laboratory information management more efficient and standardized.

Comprehensive water quality monitoring and alarm

- 3 water quality sensors, to monitor water quality and alarm (Feed water, RO water, DI water or UP water^[1]), electrode constant-0.01cm-1, temperature sensitivity-0.1 °C, and the conductivity/resistivity after temperature compensation and water temperature can be displayed simultaneously.
- 2 flow sensors, to achieve (RO water, DI water or UP water^[2]) quantitative dispensing function.
- With real-time display function of ion rejection rate of RO membrane.
- Optional feed water pressure sensor, to monitor and display feed pressure real-timely.
- Optional pure water tank liquid level sensor with dynamic icon on the main interface, can display the tank liquid level, storage capacity (L) and storage percentage (%) in real time. A glance for storage status.
- Optional independent TOC detection module, to monitor TOC of ultrapure water real-timely, detection range: 0.5-999.9ppb, detection accuracy: ±0.1ppb, in line with USP and EP system adaptability test.

^[1] According to different model, grade and type of pure water are different. For details, refer to the product manual.

^[2] According to different model, pure water type of quantitative dispensing is different. For details, refer to the product manual.



Flexible and diverse water dispensing mode^[1]

- The host is equipped with RO, DI or UP^[1] 2 standard pure water outlet. 2 kinds of water dispensing - general and quantitative spreading in parallel, to make dispensing more convenient.
- Pressure water tank, which can isolate air, stores RO water to cope with large water dispensing needs.
- Optional professional PE pure water tank^[2], to provide the third pure water outlet, to improve water dispensing efficiency.
- Up to 5 water dispenser arms per host is optional, with general, quantitative and instant - 3 kinds of water dispensing mode. More flexible to dispense.

^[1] According to different model, pure water kind and dispensing mode are different. For details, refer to the product manual.

^[2] Standard configuration of FEUE/FEDE includes 100 liters professional PE pure water tank.

Full security protection

- With DC24V as the main power supply, fully use weak current components, to meet the safety standards.
- 2 level permission management, administrator and ordinary users have strict permission distinction.
- With alarm protection of No feed water, low feed water pressure, system high pressure and full tank.
- With alarm function of standard-exceeding of Feed/RO/DI/UP water quality, and end of cartridge life.
- All alarm information can be stored in the host and cloud platform, to meet data security requirements.
- Timed shutdown and standby sleep function, to provide more protection for unattended operation.
- Optional external water leakage protection alarm device, to avoid the risk of water leakage and provide more safety protection.

Perfect consumables management

- Function of consumable life management (PP/PC/RO/DI/EDI/UP/UV/UF/TF^[1]) combined with water quality, time and capacity, to reduce consumable costs.
- The remaining life of cartridges will display on the main interface with the form of dynamic icon. Function of cartridge life reminder, step by step with 3 kinds of color, blue (normal), yellow (warning), red (alarm), all cartridges status is clearly visible.
- Click the cartridge icon to learn remaining life of cartridges and view the cartridge status timely.
- Serial number verification function of original cartridges, encrypted long serial number verification code, to prevent misoperation of cartridges installation and replacement.


^[1] According to different model, cartridges are different. For details, refer to the product manual.

Advanced System Configuration Guaranteeing Strong Quality



1 Powerful 20-inch pretreated cartridge

- PP cartridge with deep folding membrane, accuracy of $5\mu\text{m}$, to filter particles efficiently in source water.
- PC cartridge with high performance activated carbon fiber with catalyst, accuracy of $5\mu\text{m}$, to adsorb organics and residual chlorine efficiently and avoid carbon powder precipitation maximumly.
- The combination of folding filter & carbon fiber filter can bring greater cartridge processing capacity, extend the replacement cycle and reduce the running cost.

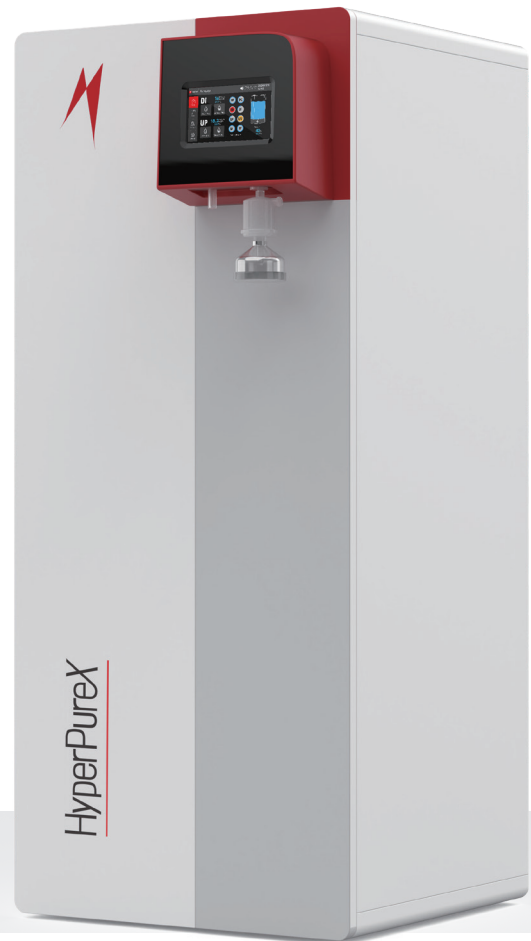

2 Rigorous double RO system^[1]

- Double RO system can remove up to 99% soluble inorganic ions, 99% soluble organics, microorganisms and particles.
- Compared with single RO system, the double RO water quality can be stable $< 5\mu\text{S/cm}$ (feed water conductivity $< 1500\mu\text{S/cm}$), and the life of the ultrapure unit is longer.
- Equipped with high quality RO membrane, to achieve combination of long life, stability and high ion rejection rate.
- Auto-flushing function of RO module with adjustable flushing interval and duration, to effectively prevent scale and prolong the life of the membrane.
- The automatic discharge function of unqualified RO water can ensure that the RO water quality is suitable to enter the back-end module.
- Integral package of discarded RO module, easy to install and maintain.



3 High performance purification cartridge

- Patented cartridge structure can prevent the stratification of resin and ensure the exchange capacity of cartridge.
- The resin filling capacity per cartridge is up to 5.2 liters, and up to 2 cartridges^[1] can be equipped every host, with a total filling capacity of 10.4 liters, achieving greater ion exchange capacity and significantly reducing the running cost.
- All DuPont resin and high purity material of column ensure absolute $18.2\text{M}\Omega\cdot\text{cm}$ of ultrapure water resistivity and reduce TOC precipitation.

^[1] When system output $\geq 150\text{ l/h}$, the purification resin cartridge uses external super large DI cartridge with capacity of up to 32L resin filling.

^[1] Applicable to FEUS/FEDS/FERS.



4 Built-in 7.5-liter pressure water tank

- With dual functions of water storage and pressurization, FDA approved, its fully enclosed structure effectively isolates air, and prevent the touching of CO₂ and other pollutants with pure water. Up to 100 liters is optional volume.
- Optional 100 / 250 liters professional PE pure water tank, equipped with liquid level sensor and air filter, to meet the needs of large water consumption.




5 Professional mini EDI module^[1]

- Without softener and chemical regeneration, pure water in grade II, with resistivity > 10MΩ.cm@25 °C (generally above 15MΩ.cm), and TOC<30ppb^[2], is available. It can prolong the life of back-end module and reduce the running cost.

^[1] Applicable to FEUE/FEDE.

^[2] The values vary depending on the nature and concentration of contaminants in source water.

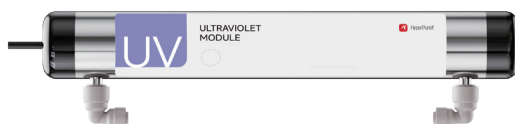

6 Double wavelength UV module^[1]

- Long-life ultraviolet lamp (185&254nm), combined with SUS316L flow shell, can reduce the value of TOC to ≤2ppb^[2], and can achieve efficient sterilization and inhibit bacterial growth, suitable for HPLC, UPLC, LC-MS and other precision instruments.

^[1] Applicable to ultrapure water systems equipped with UV module.

^[2] The values vary depending on the nature and concentration of contaminants in source water.



7 Ultrafiltration module^[1]

- With PES membrane and MWCO>5000D, effectively removes pyrogen/endotoxin, RNase, DNase, and produces nuclease-free, proteinase-free and bacterial-free ultrapure water, suitable for life science applications, such as cell culture/IVF.

^[1] Applicable to ultrapure water systems equipped with UF module.



8 MF terminal microfilter^[1]

- (0.45±0.2)μm double-layer PES membrane ensures microbial retention, effectively removes particles and bacteria, and meets critical application requirements.

^[1] Applicable to ultrapure or high-pure water system. For details, refer to the product manual.



9 UF terminal ultrafilter^[1]

- With PES membrane and MWCO>15000D, effectively removes pyrogen/endotoxin, RNase, DNase, and produces nuclease-free, proteinase-free and bacterial-free ultrapure water, suitable for life science applications, such as cell culture/IVF.

^[1] Optional accessory for ultrapure water system only.



10 Multiple communication interfaces^[1]

- USB interface, to export running data or upgrade system version online.
- RJ45 or USB/WIFI interface, to achieve the IOT and cloud platform connection.
- HiDis pure water dispensing arm interface, to achieve power supply and data communication with the host. (Optional)
- L-Tank pure water tank interface, to synchronize water tank level signal with the host. (Optional)
- FS foot switch interface, easy to dispense, suitable for more dispensing scenarios. (Optional)

^[1] The interface configuration varies according to the model. For details, refer to the product manual.

FEU Eminent series

Medium Ultrapure Water System


— Ultrapure water, high pure water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, stable and reliable single RO system, and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank.

System output: 60, 120, 180, 250 liters/h. It can simultaneously produce ultrapure water (18.2MΩ.cm) and high pure water (>16MΩ.cm). The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ASTM D1193-06, GB/T 11446.1-2013, GB/T 33087-2016, GB/T 6682-2008, CP, EP, USP, JP, CAP, CLSI, etc.



Flow Diagram

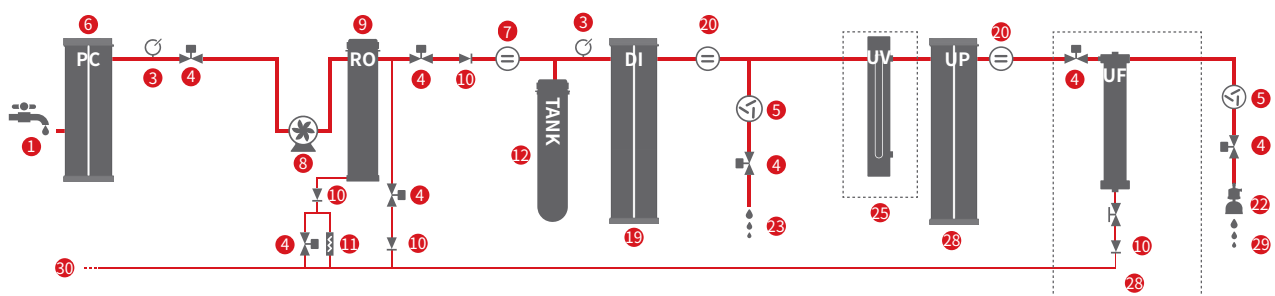

- | | | | |
|-----------------------------|-----------------------|-----------------------|--------------------------|
| ① Feed Water | ⑨ RO cartridge | ⑰ Three way valve | ⑳ UV Component |
| ② PP Pretreatment Cartridge | ⑩ One way valve | ⑱ High tension switch | ㉑ TOC Component |
| ③ Pressure sensor | ⑪ Flow Restrictor | ㉒ DI Cartridge | ㉒ UP Ultrapure cartridge |
| ④ Solenoid valve | ⑫ Pressure water tank | ㉓ Resistivity Sensor | ㉓ UF Cartridge |
| ⑤ Flow sensor | ⑬ RO Water Outlet | ㉔ Sanitization Block | ㉔ UP Water Outlet |
| ⑥ PC Pretreatment Cartridge | ⑭ Low tension switch | ㉕ Final Filter | ㉕ Drain Outlet |
| ⑦ Conductivity Sensor | ⑮ EDI Component | ㉖ DI Water Outlet | |
| ⑧ Pump | ⑯ PE water tank | ㉖ Dispenser arm | |

FEU Specifications

Name	Standard	Low TOC	Eliminating endotoxin	Synthesizing
Model	FEU-60/120/180/250	FEU-60/120/180/250UV	FEU-60/120/180/250UF	FEU-60/120/180/250UVF
Production rate ^[1]	60 series: 60 L/hour, 120 series: 120 L/hour, 180 series: 180 L/hour, 250 series: 250 L/hour			
Dispensing rate ^[2]	Up to 4.2 liters/minute	Up to 4.2 liters/minute	Up to 4.2 liters/minute	Up to 4.2 liters/minute
Ultrapure water quality ^[3]				
Resistivity (25°C) ^[4]	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm
Conductivity (25°C)	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm
TOC ^[5]	10 ppb ^[6]	3 ppb ^[7]	10 ppb ^[6]	3 ppb ^[7]
Particles ^[8]	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)
Bacteria ^[9]	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml
Endotoxin ^[10]	N/A	N/A	<0.001 EU/ml	<0.001 EU/ml
RNases ^[10]	N/A	N/A	1 pg/ml	1 pg/ml
DNases ^[10]	N/A	N/A	5 pg/ml	5 pg/ml
Protease ^[10]	N/A	N/A	0.15 μg/ml	0.15 μg/ml
DI water quality ^[3]				
Resistivity (25°C) ^[4]	>17.5 MΩ.cm	>17.5 MΩ.cm	>17.5 MΩ.cm	>17.5 MΩ.cm
Conductivity (25°C)	<0.057 μs/cm	<0.057 μs/cm	<0.057 μs/cm	<0.057 μs/cm
Particles ^[8]	N/A	N/A	N/A	N/A
Bacteria ^[9]	N/A	N/A	N/A	N/A
RO ^{1st} water quality ^[3]				
Ion rejection rate	98%-99% (with new RO module)	98%-99% (with new RO module)	98%-99% (with new RO module)	98%-99% (with new RO module)
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%	>99%	>99%
Feed water requirements				
Water source type	Tap water	Tap water	Tap water	Tap water
Pressure	1-6 bar	1-6 bar	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm	<3 ppm	<3 ppm
PH	4-10	4-10	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz
Total Power	60series: 240W; 120series: 300W; 180series: 400W; 250series: 600W			
Dimension (L×W×H)	450×521×1016mm	450×521×1016mm	450×521×1016mm	450×521×1016mm
Weight	About 64kg	About 64kg	About 64kg	About 64kg
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FED Eminent series

Medium Pure Water System


— High pure water, RO^{1st} water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, stable and reliable single RO system, and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank.

System output: 60, 120, 180, 250 liters/h. It can simultaneously produce high pure water (>17.5MΩ.cm) and single RO water. The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ISO3696 (Grade 2), GB/T 6682 (Grade 1), ASTM D1193 (Type II reagent water), JIS K0557, etc., also meets the purified water technical requirements of CP, EP, USP, JP and other national pharmacopoeia.



Flow Diagram

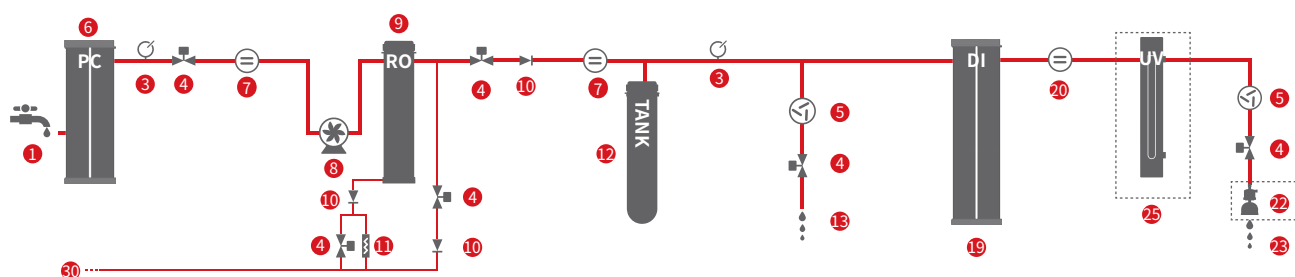

- | | | | |
|-----------------------------|-----------------------|-----------------------|--------------------------|
| ① Feed Water | ⑨ RO cartridge | ⑰ Three way valve | ⑵ UV Component |
| ② PP Pretreatment Cartridge | ⑩ One way valve | ⑱ High tension switch | ⑶ TOC Component |
| ③ Pressure sensor | ⑪ Flow Restrictor | ⑲ DI Cartridge | ⑴ UP Ultrapure cartridge |
| ④ Solenoid valve | ⑫ Pressure water tank | ⑳ Resistivity Sensor | ⑵ UF Cartridge |
| ⑤ Flow sensor | ⑬ RO Water Outlet | ㉑ Sanitization Block | ⑶ UP Water Outlet |
| ⑥ PC Pretreatment Cartridge | ⑭ Low tension switch | ㉒ Final Filter | ⑷ Drain Outlet |
| ⑦ Conductivity Sensor | ⑮ EDI Component | ㉓ DI Water Outlet | |
| ⑧ Pump | ⑯ PE water tank | ㉔ Dispenser arm | |

FED Specifications

Name	Standard	Eliminating bacteria and particle
Model	FED-60/120/180/250	FED-60/120/180/250UT
Production rate ^[1]	60 series: 60 L/hour, 120 series: 120 L/hour, 180 series: 180 L/hour, 250 series: 250 L/hour	
Dispensing rate ^[2]	Up to 4.2 liters/minute	Up to 4.2 liters/minute
DI water quality ^[3]		
Resistivity (25°C) ^[4]	>17.5 MΩ.cm	>17.5 MΩ.cm
Conductivity (25°C)	<0.057 μs/cm	<0.057 μs/cm
Particles ^[8]	N/A	<1 /ml (>0.2μm)
Bacteria ^[9]	N/A	<0.01 CFU/ml
RO^{1st} water quality ^[3]		
Ion rejection rate	98%-99% (with new RO module)	98%-99% (with new RO module)
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%
Feed water requirements		
Water source type	Tap water	Tap water
Pressure	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO₃)	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm
PH	4-10	4-10
Dissolved CO₂	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz
Total Power	60series:240W, 120series:300W 180series:400W, 250series:600W	60series:240W, 120series:300W 180series:400W, 250series:600W
Dimension (L×W×H)	450×521×1016mm	450×521×1016mm
Weight	About 63KG	About 63KG
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FEUS Eminent series

Medium Ultrapure Water System


— Ultrapure water, high pure water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, rigorous double RO system, and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank.

System output: 50, 100, 150 liters/h. It can simultaneously produce ultrapure water (18.2MΩ.cm) and high pure water (>16MΩ.cm). The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ASTM D1193-06, GB/T 11446.1-2013, GB/T 33087-2016, GB/T 6682-2008, CP, EP, USP, JP, CAP, CLSI, etc.



Flow Diagram

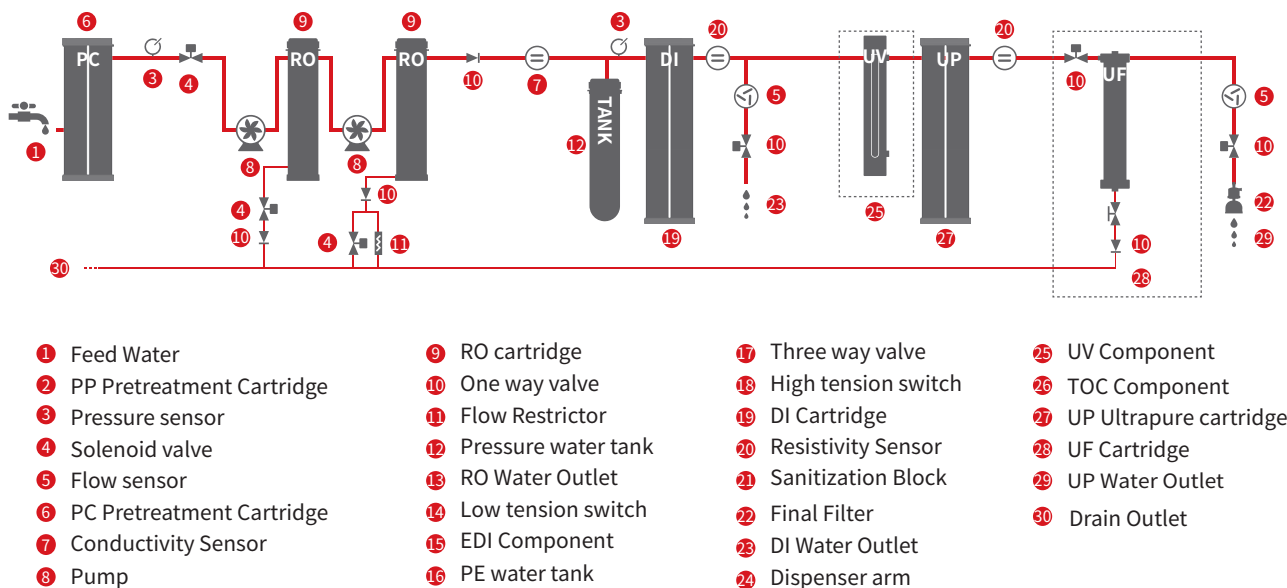

FEUS Specifications

Name	Standard	Low TOC	Eliminating endotoxin	Synthesizing
Model	FEUS-50/100/150	FEUS-50/100/150UV	FEUS-50/100/150UF	FEUS-50/100/150UVF
Production rate ^[1]	50series: 50 L/hour; 100series: 100 L/hour; 150series: 150 L/hour			
Dispensing rate ^[2]	Up to 2.5 liters/minute	Up to 2.5 liters/minute	Up to 2.5 liters/minute	Up to 2.5 liters/minute
Ultrapure water quality ^[3]				
Resistivity (25°C) ^[4]	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm
Conductivity (25°C)	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm
TOC ^[5]	10 ppb ^[6]	3 ppb ^[7]	10 ppb ^[6]	3 ppb ^[7]
Particles ^[8]	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)
Bacteria ^[9]	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml
Endotoxin ^[10]	N/A	N/A	<0.001 EU/ml	<0.001 EU/ml
RNases ^[10]	N/A	N/A	1 pg/ml	1 pg/ml
DNases ^[10]	N/A	N/A	5 pg/ml	5 pg/ml
Protease ^[10]	N/A	N/A	0.15 μg/ml	0.15 μg/ml
DI water quality ^[3]				
Resistivity (25°C) ^[4]	>17.5 MΩ.cm	>17.5 MΩ.cm	>17.5 MΩ.cm	>17.5 MΩ.cm
Conductivity (25°C)	<0.057 μs/cm	<0.057 μs/cm	<0.057 μs/cm	<0.057 μs/cm
Particles ^[8]	N/A	N/A	N/A	N/A
Bacteria ^[9]	N/A	N/A	N/A	N/A
RO ^{2nd} water quality ^[3]				
Resistivity (25°C) ^[4]	>0.2 MΩ.cm	>0.2 MΩ.cm	>0.2 MΩ.cm	>0.2 MΩ.cm
Conductivity (25°C)	<5 μs/cm	<5 μs/cm	<5 μs/cm	<5 μs/cm
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%	>99%	>99%
Feed water requirements				
Water source type	Tap water	Tap water	Tap water	Tap water
Pressure	1-6 bar	1-6 bar	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm	<3 ppm	<3 ppm
PH	4-10	4-10	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz
Total Power	50series: 240W; 100series: 400W; 150series: 600W			
Dimension (L×W×H)	450×521×1016mm	450×521×1016mm	450×521×1016mm	450×521×1016mm
Weight	About 70KG	About 70KG	About 70KG	About 70KG
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FEDS Eminent series

Medium Pure Water System


—— High pure water, RO^{2nd} water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, rigorous double RO system, and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank.

System output: 50, 100, 150 liters/h. It can simultaneously produce high pure water (>17.5MΩ.cm) and double RO water (<5μs/cm). The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ISO3696 (Grade 2), GB/T 6682 (Grade 1), ASTM D1193 (Type II reagent water), JIS K0557, etc., also meets the purified water technical requirements of CP, EP, USP, JP and other national pharmacopoeia.



Flow Diagram

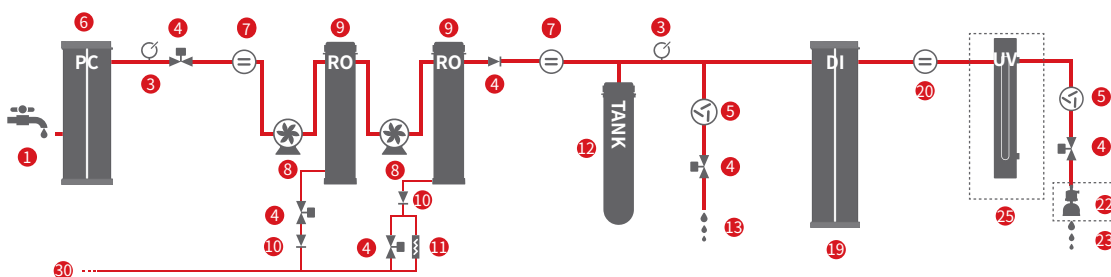

- | | | | |
|-----------------------------|-----------------------|-----------------------|--------------------------|
| ① Feed Water | ⑨ RO cartridge | ⑰ Three way valve | ⑳ UV Component |
| ② PP Pretreatment Cartridge | ⑩ One way valve | ⑱ High tension switch | ㉑ TOC Component |
| ③ Pressure sensor | ⑪ Flow Restrictor | ㉒ DI Cartridge | ㉒ UP Ultrapure cartridge |
| ④ Solenoid valve | ⑫ Pressure water tank | ㉓ Resistivity Sensor | ㉓ UF Cartridge |
| ⑤ Flow sensor | ⑬ RO Water Outlet | ㉔ Sanitization Block | ㉔ UP Water Outlet |
| ⑥ PC Pretreatment Cartridge | ⑭ Low tension switch | ㉕ Final Filter | ㉕ Drain Outlet |
| ⑦ Conductivity Sensor | ⑮ EDI Component | ㉖ DI Water Outlet | |
| ⑧ Pump | ⑯ PE water tank | ㉖ Dispenser arm | |

FEDS Specifications

Name	Standard	Eliminating bacteria and particle
Model	FEDS-50/100/150	FEDS-50/100/150UT
Production rate ^[1]	50series: 50 L/hour; 100series: 100 L/hour; 150series: 150 L/hour	
Dispensing rate ^[2]	Up to 2.5 liters/minute	Up to 2.5 liters/minute
DI water quality ^[3]		
Resistivity (25°C) ^[4]	>17.5 MΩ.cm	>17.5 MΩ.cm
Conductivity (25°C)	<0.057 μs/cm	<0.057 μs/cm
Particles ^[8]	N/A	<1 /ml (>0.2μm)
Bacteria ^[9]	N/A	<0.01 CFU/ml
RO ^{2nd} water quality ^[3]		
Resistivity (25°C) ^[4]	>0.2 MΩ.cm	>0.2 MΩ.cm
Conductivity (25°C)	<5 μs/cm	<5 μs/cm
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%
Feed water requirements		
Water source type	Tap water	Tap water
Pressure	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm
PH	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz
Total Power	50series: 240W; 100series: 400W; 150series: 600W	50series: 240W; 100series: 400W; 150series: 600W
Dimension (L×W×H)	450×521×1016mm	450×521×1016mm
Weight	About 69KG	About 69KG
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FEUE Eminent series

Medium Ultrapure Water System


— Ultrapure water, EDI water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, stable and rigorous double RO system, advanced EDI module and DI ion-exchange cartridges with larger capacity, equipping with built-in 7.5-liter pressure water tank and professional-grade pure water tank with 100-liter.

System output: 50, 100, 150 liters/h. Maximum output per day is up to 3600 liters. It can simultaneously produce ultrapure water ($18.2\text{M}\Omega\cdot\text{cm}$) and EDI water (Resistivity $>10\text{M}\Omega\cdot\text{cm}$, $\text{TOC} < 30\text{ppb}$) with optimized running cost. The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ASTM D1193-06, GB/T 11446.1-2013, GB/T 33087-2016, GB/T 6682-2008, CP, EP, USP, JP, CAP, CLSI, etc.



Flow Diagram

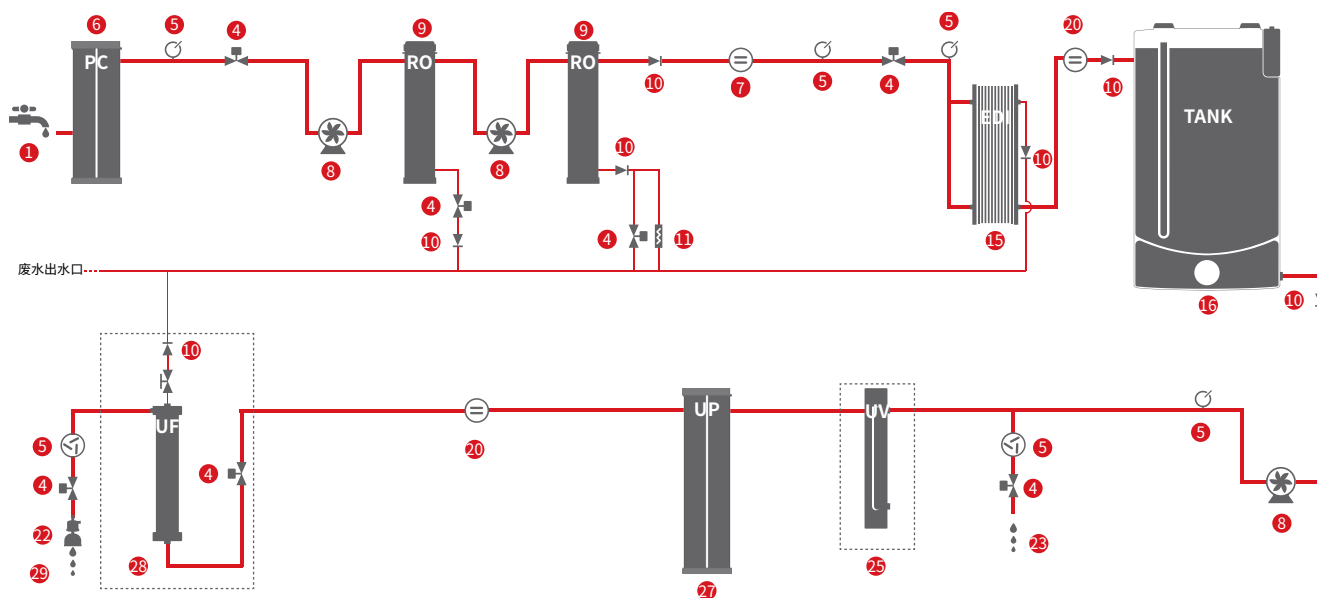

- | | | | |
|-----------------------------|-----------------------|-----------------------|---------------------------|
| ① Feed Water | ⑨ RO cartridge | ⑰ Three way valve | ②⑤ UV Component |
| ② PP Pretreatment Cartridge | ⑩ One way valve | ⑱ High tension switch | ②⑥ TOC Component |
| ③ Pressure sensor | ⑪ Flow Restrictor | ⑲ DI Cartridge | ②⑦ UP Ultrapure cartridge |
| ④ Solenoid valve | ⑫ Pressure water tank | ⑳ Resistivity Sensor | ②⑧ UF Cartridge |
| ⑤ Flow sensor | ⑬ RO Water Outlet | ㉑ Sanitization Block | ②⑨ UP Water Outlet |
| ⑥ PC Pretreatment Cartridge | ⑭ Low tension switch | ㉒ Final Filter | ③⑩ Drain Outlet |
| ⑦ Conductivity Sensor | ⑮ EDI Component | ㉓ DI Water Outlet | |
| ⑧ Pump | ⑯ PE water tank | ㉔ Dispenser arm | |

FEUE Specifications

Name	Standard	Low TOC	Eliminating endotoxin	Synthesizing
Model	FEUE-50/100/150	FEUE-50/100/150UV	FEUE-50/100/150UF	FEUE-50/100/150UVF
Production rate ^[1]	50series: 50 L/hour; 100series: 100 L/hour; 150series: 150 L/hour			
Dispensing rate ^[2]	Up to 2.5 liters/minute	Up to 2.5 liters/minute	Up to 2.5 liters/minute	Up to 2.5 liters/minute
Ultrapure water quality ^[3]				
Resistivity (25°C) ^[4]	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm	18.2 MΩ.cm
Conductivity (25°C)	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm	0.055 μs/cm
TOC ^[5]	10 ppb ^[6]	3 ppb ^[7]	10 ppb ^[6]	3 ppb ^[7]
Particles ^[8]	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)	<1 /ml (>0.2μm)
Bacteria ^[9]	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml	<0.01 CFU/ml
Endotoxin ^[10]	N/A	N/A	<0.001 EU/ml	<0.001 EU/ml
RNases ^[10]	N/A	N/A	1 pg/ml	1 pg/ml
DNases ^[10]	N/A	N/A	5 pg/ml	5 pg/ml
Protease ^[10]	N/A	N/A	0.15 μg/ml	0.15 μg/ml
EDI water quality ^[3]				
Resistivity (25°C) ^[4]	>10 MΩ.cm	>10 MΩ.cm	>10 MΩ.cm	>10 MΩ.cm
Conductivity (25°C)	<0.1 μs/cm	<0.1 μs/cm	<0.1 μs/cm	<0.1 μs/cm
TOC ^[5]	≤ 30 ppb	≤ 30 ppb	≤ 30 ppb	≤ 30 ppb
Particles ^[8]	N/A	N/A	N/A	N/A
Bacteria ^[9]	N/A	N/A	N/A	N/A
RO ^{2nd} water quality ^[3]				
Resistivity (25°C) ^[4]	>0.2 MΩ.cm	>0.2 MΩ.cm	>0.2 MΩ.cm	>0.2 MΩ.cm
Conductivity (25°C)	<5 μs/cm	<5 μs/cm	<5 μs/cm	<5 μs/cm
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%	>99%	>99%
Feed water requirements				
Water source type	Tap water	Tap water	Tap water	Tap water
Pressure	1-6 bar	1-6 bar	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm	<3 ppm	<3 ppm
PH	4-10	4-10	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz	220V, 50/60Hz
Total Power	50series: 240W; 100series: 400W; 150series: 600W			
Dimension (L×W×H)	Main host: 450×521×1016mm Tank: 450×300×1016mm	Main host: 450×521×1016mm Tank: 450×300×1016mm	Main host: 450×521×1016mm Tank: 450×300×1016mm	Main host: 450×521×1016mm Tank: 450×300×1016mm
Weight	About 75KG	About 75KG	About 75KG	About 75KG
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FEDE Eminent series

Medium Pure Water System


——EDI water, RO^{2nd} water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure, stable and rigorous double RO system, and advanced EDI module, equipping with built-in 7.5-liter pressure water tank and professional-grade pure water tank with 100-liter.

System output: 50, 100, 150 liters/h. Maximum output per day is up to 3600 liters. It can simultaneously produce double RO water ($<5\mu\text{s/cm}$) and EDI water (Resistivity $>10\text{M}\Omega\cdot\text{cm}$, TOC $<30\text{ppb}$) with optimized running cost. The quality of pure water fully meets or exceeds the requirements of water quality standard specified by ISO3696 (Grade 2), GB/T 6682 (Grade 1), ASTM D1193 (Type II reagent water), JIS K0557, etc., also meets the purified water technical requirements of CP, EP, USP, JP and other national pharmacopoeia.



Flow Diagram

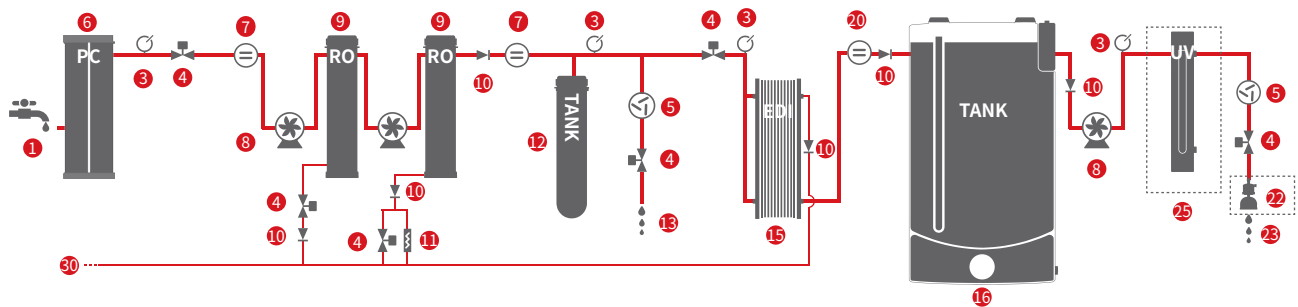

- | | | | |
|-----------------------------|-----------------------|-----------------------|---------------------------|
| ① Feed Water | ⑨ RO cartridge | ⑰ Three way valve | ②⑤ UV Component |
| ② PP Pretreatment Cartridge | ⑩ One way valve | ⑱ High tension switch | ②⑥ TOC Component |
| ③ Pressure sensor | ⑪ Flow Restrictor | ⑲ DI Cartridge | ②⑦ UP Ultrapure cartridge |
| ④ Solenoid valve | ⑫ Pressure water tank | ⑳ Resistivity Sensor | ②⑧ UF Cartridge |
| ⑤ Flow sensor | ⑬ RO Water Outlet | ㉑ Sanitization Block | ②⑨ UP Water Outlet |
| ⑥ PC Pretreatment Cartridge | ⑭ Low tension switch | ㉒ Final Filter | ③⑩ Drain Outlet |
| ⑦ Conductivity Sensor | ⑮ EDI Component | ㉓ DI Water Outlet | |
| ⑧ Pump | ⑯ PE water tank | ㉔ Dispenser arm | |

FEDE Specifications

Name	Standard	Eliminating bacteria and particle
Model	FEDE-50/100/150	FEDE-50/100/150UT
Production rate ^[1]	50series: 50 L/hour; 100series: 100 L/hour; 150series: 150 L/hour	
Dispensing rate ^[2]	Up to 2.5 liters/minute	Up to 2.5 liters/minute
EDI water quality ^[3]		
Resistivity (25°C) ^[4]	>10 MΩ.cm	>10 MΩ.cm
Conductivity (25°C)	<0.1 μs/cm	<0.1 μs/cm
TOC ^[5]	≤ 30 ppb	≤ 30 ppb
Particles ^[8]	N/A	<1 /ml (>0.2μm)
Bacteria ^[9]	N/A	<0.01 CFU/ml
RO ^{2nd} water quality ^[3]		
Resistivity (25°C) ^[4]	>0.2 MΩ.cm	>0.2 MΩ.cm
Conductivity (25°C)	<5 μs/cm	<5 μs/cm
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%
Feed water requirements		
Water source type	Tap water	Tap water
Pressure	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm
PH	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz
Total Power	50series: 240W; 100series: 400W; 150series: 600W	50series: 240W; 100series: 400W; 150series: 600W
Dimension (L×W×H)	Main host: 450×521×1016mm Tank: 450×300×1016mm	Main host: 450×521×1016mm Tank: 450×300×1016mm
Weight	About 68kg	About 68kg
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set 100-liter PE water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

FERS Eminent series

Medium Double RO Water System


——RO^{2nd} water, RO^{1st} water

With tap water inlet, using the innovative human-computer interactive control system and 5-inch IPS colorful capacitive touch screen, integrating functions of Internet of Things (IOT) and cloud platform, embedding new purification cartridges with patented structure and rigorous double RO system, equipping with built-in 1.8-liter pressure water tank.

System output: 50, 100, 150 liters/h. It can simultaneously produce single RO and double RO water. The ion rejection rate of single RO water is above of 98%, and the conductivity of double RO water is less than 5 μ s/cm. The quality of pure water fully meets or exceeds the requirements of water quality standard specified by GB/T 6682-2008 (Grade 3).



Flow Diagram

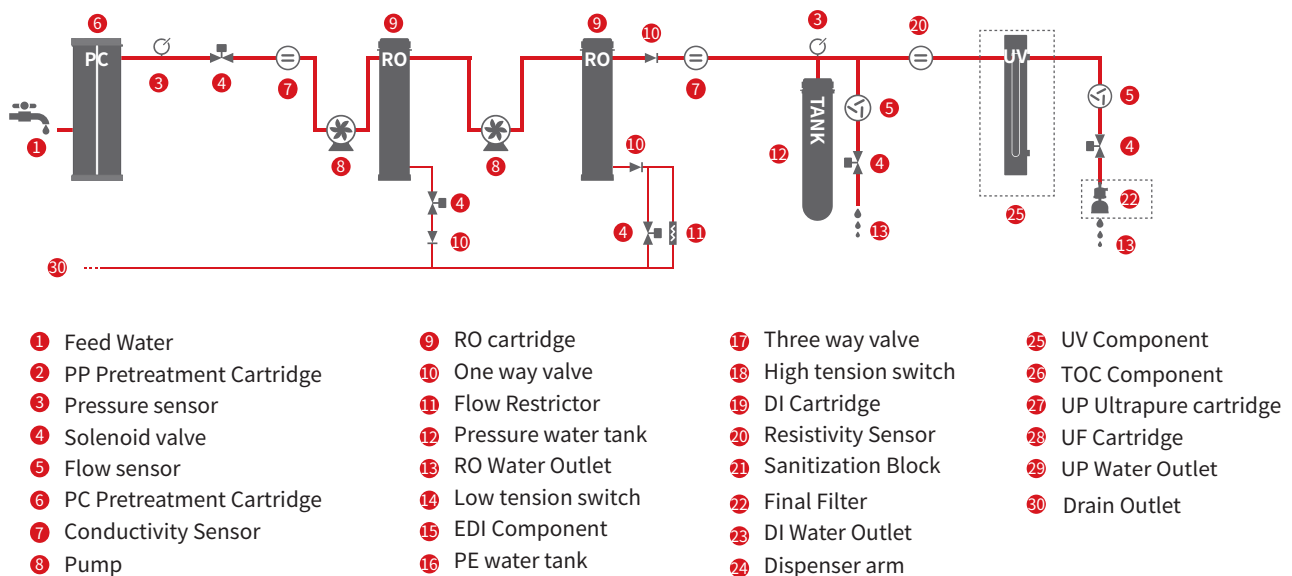

FERS Specifications

Name	Standard	Eliminating bacteria and particle
Model	FERS-50/100/150	FERS-50/100/150UT
Production rate ^[1]	50series: 50 L/hour; 100series: 100 L/hour; 150series: 150 L/hour	
Dispensing rate ^[2]	Up to 2.5 liters/minute	Up to 2.5 liters/minute
RO ^{1st} water quality ^[3]		
Ion rejection rate	>98% (with new RO module)	>98% (with new RO module)
RO ^{2nd} water quality ^[3]		
Resistivity (25°C) ^[4]	>0.2 MΩ.cm	>0.2 MΩ.cm
Conductivity (25°C)	<5 μs/cm	<5 μs/cm
Organic rejection rate	>99% (MW>300 Dalton)	>99% (MW>300 Dalton)
Particles and bacteria rejection rate	>99%	>99%
Particles ^[8]	N/A	<1 /ml (>0.2μm)
Bacteria ^[9]	N/A	<0.01 CFU/ml
Feed water requirements		
Water source type	Tap water	Tap water
Pressure	1-6 bar	1-6 bar
Temperature	5-40 °C	5-40 °C
Conductivity	<2000 μs/cm	<2000 μs/cm
Total hardness (In CaCO ₃)	<300 ppm	<300 ppm
TOC	<2000 ppb	<2000 ppb
Free chlorine	<3 ppm	<3 ppm
PH	4-10	4-10
Dissolved CO ₂	<40 ppm	<40 ppm
Power supply	220V, 50/60Hz	220V, 50/60Hz
Total Power	50series: 240W; 100series: 400W; 150series: 600W	50series: 240W; 100series: 400W; 150series: 600W
Dimension (L×W×H)	450×521×1016mm	450×521×1016mm
Weight	About 60kg	About 60kg
Standard configuration	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set	Main host 1 set All cartridges 1 set 7.5-liter pressure water tank 1 set

[1] Affected by inlet water quality, pressure, temperature and status of RO membrane

[2] Affected by the tank status and terminal filter

[3] The following values are typical and may vary depending on the nature and concentration of feed water contaminants

[4] According to USP requirements, the resistivity can be displayed as a non-temperature-compensated value

[5] Affected by the type of organics

[6] Inlet TOC<1000ppb, follow professional operating procedures and correct sampling conditions

[7] Inlet TOC<50ppb, follow professional operating procedures and correct sampling conditions

[8] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[9] Equip with terminal microfilter and follow professional operating procedures and correct sampling conditions

[10] Equip with terminal ultrafilter and follow professional operating procedures and correct sampling conditions

Hyperpurex® lab water system

Bring you products and services beyond expected

ISO
3696
US Pharmacopoeia
GB/T 33087 2016
Japan Pharmacopoeia
ISO9001
CLSI GB/T.11446 1-2013
ASTM D 5196 **ISO14001**
China Pharmacopoeia ASTM
GB/T 6682-2008 JIS K 0557
Eu Pharmacopoeia D1193
CE Quality Standard
HyperpureX®

PRODUCT

- Under management system of ISO9001 and ISO14001, in accordance with CE quality standards, we carry out product design, research & development and manufacturing to ensure long-term stability and reliability of quality.
- To help you meet industry specifications, we can assist in providing certificates of conformity, calibration certificates, quality certificates, performance reports, water quality compliance certificates and other supporting documents upon request.
- Hyperpurex® FE Eminent series - lab water system can produce pure water/ultrapure water to meet the requirements of the following organizations:
- Chinese Pharmacopoeia-CP, United States Pharmacopoeia-USP, European Pharmacopoeia-EP, Japanese Pharmacopoeia -JP, GB/T 33087-2016, GB/T 6682-2008, GB/T 11446.1-2013, ASTM D1193, ASTM D 5196, ISO 3696, CLSI, JIS K 0557.

SERVICE

We wholeheartedly serve, only for your full satisfaction.

With customer satisfaction as the service goal, to continue to create value for customers as the direction, to grow together with customers as the concept, based on professionalism, we are full of sincerity and enthusiasm, committing to providing customers with professional and perfect technical support and after-sales service. So that you can devote all your energy to focus on the work.

Our service include:

- 24 months product warranty (excluding filter consumables)
- On-site professional training of installation, use and maintenance.
- Regular engineer return visit service
- Free continuous optimization and upgrading service of product life cycle.
- Professional and rigorous 3Q(IQ/OQ/PQ) verification documentation and verification services in both English and Chinese, to help you meet compliance requirements of GLP, GMP and cGMP.

Ordering Information

Host	Item No	Product description
	FEU-60	Medium ultrapure water system,60L/h, Standard, Ultrapure water, high pure water
	FEU-120	Medium ultrapure water system,120L/h, Standard, Ultrapure water, high pure water
	FEU-180	Medium ultrapure water system,180L/h, Standard, Ultrapure water, high pure water
	FEU-250	Medium ultrapure water system,250L/h, Standard, Ultrapure water, high pure water
	FEU-60UV	Medium ultrapure water system,60L/h, Low TOC, Ultrapure water, high pure water
	FEU-120UV	Medium ultrapure water system,120L/h, Low TOC, Ultrapure water, high pure water
	FEU-180UV	Medium ultrapure water system,180L/h, Low TOC, Ultrapure water, high pure water
	FEU-250UV	Medium ultrapure water system,250L/h, Low TOC, Ultrapure water, high pure water
	FEU-60UF	Medium ultrapure water system,60L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEU-120UF	Medium ultrapure water system,120L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEU-180UF	Medium ultrapure water system,180L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEU-250UF	Medium ultrapure water system,250L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEU-60UVF	Medium ultrapure water system,60L/h, Synthesizing, Ultrapure water, high pure water
	FEU-120UVF	Medium ultrapure water system,120L/h, Synthesizing, Ultrapure water, high pure water
	FEU-180UVF	Medium ultrapure water system,180L/h, Synthesizing, Ultrapure water, high pure water
	FEU-250UVF	Medium ultrapure water system,250L/h, Synthesizing, Ultrapure water, high pure water
	FED-60	Medium pure water system,60L/h, Standard, High pure water, RO ^{1st} water
	FED-120	Medium pure water system,120L/h, Standard, High pure water, RO ^{1st} water
	FED-180	Medium pure water system,180L/h, Standard, High pure water, RO ^{1st} water
	FED-250	Medium pure water system,250L/h, Standard, High pure water, RO ^{1st} water
	FED-60UT	Medium pure water system,60L/h, Eliminating bacteria and particle, High pure water, RO ^{1st} water
	FED-120UT	Medium pure water system,120L/h, Eliminating bacteria and particle, High pure water, RO ^{1st} water
	FED-180UT	Medium pure water system,180L/h, Eliminating bacteria and particle, High pure water, RO ^{1st} water
	FED-250UT	Medium pure water system,250L/h, Eliminating bacteria and particle, High pure water, RO ^{1st} water
	FEUS-50	Medium ultrapure water system,50L/h, Standard, Ultrapure water, high pure water
	FEUS-100	Medium ultrapure water system,100L/h, Standard, Ultrapure water, high pure water
	FEUS-150	Medium ultrapure water system,150L/h, Standard, Ultrapure water, high pure water
	FEUS-50UV	Medium ultrapure water system,50L/h, Low TOC, Ultrapure water, high pure water
	FEUS-100UV	Medium ultrapure water system,100L/h, Low TOC, Ultrapure water, high pure water
	FEUS-150UV	Medium ultrapure water system,150L/h, Low TOC, Ultrapure water, high pure water
	FEUS-50UF	Medium ultrapure water system,50L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEUS-100UF	Medium ultrapure water system,100L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEUS-150UF	Medium ultrapure water system,150L/h, Eliminating endotoxin, Ultrapure water, high pure water
	FEUS-50UVF	Medium ultrapure water system,50L/h, Synthesizing, Ultrapure water, high pure water
	FEUS-100UVF	Medium ultrapure water system,100L/h, Synthesizing, Ultrapure water, high pure water
	FEUS-150UVF	Medium ultrapure water system,150L/h, Synthesizing, Ultrapure water, high pure water
	FEDS-50	Medium pure water system,50L/h, Standard, High pure water, RO ^{2nd} water
	FEDS-100	Medium pure water system,100L/h, Standard, High pure water, RO ^{2nd} water
	FEDS-150	Medium pure water system,150L/h, Standard, High pure water, RO ^{2nd} water
	FEDS-50UT	Medium pure water system,50L/h, Eliminating bacteria and particle, High pure water, RO ^{2nd} water
	FEDS-100UT	Medium pure water system,100L/h, Eliminating bacteria and particle, High pure water, RO ^{2nd} water
	FEDS-150UT	Medium pure water system,150L/h, Eliminating bacteria and particle, High pure water, RO ^{2nd} water
	FEUE-50	Medium ultrapure water system,50L/h, Standard, Ultrapure water, EDI water
	FEUE-100	Medium ultrapure water system,100L/h, Standard, Ultrapure water, EDI water
	FEUE-150	Medium ultrapure water system,150L/h, Standard, Ultrapure water, EDI water
	FEUE-50UV	Medium ultrapure water system,50L/h, Low TOC, Ultrapure water, EDI water
	FEUE-100UV	Medium ultrapure water system,100L/h, Low TOC, Ultrapure water, EDI water
	FEUE-150UV	Medium ultrapure water system,150L/h, Low TOC, Ultrapure water, EDI water
	FEUE-50UF	Medium ultrapure water system,50L/h, Eliminating endotoxin, Ultrapure water, EDI water
	FEUE-100UF	Medium ultrapure water system,100L/h, Eliminating endotoxin, Ultrapure water, EDI water
	FEUE-150UF	Medium ultrapure water system,150L/h, Eliminating endotoxin, Ultrapure water, EDI water
	FEUE-50UVF	Medium ultrapure water system,50L/h, Synthesizing, Ultrapure water, EDI water
	FEUE-100UVF	Medium ultrapure water system,100L/h, Synthesizing, Ultrapure water, EDI water
	FEUE-150UVF	Medium ultrapure water system,150L/h, Synthesizing, Ultrapure water, EDI water
	FEDE-50	Medium pure water system,50L/h, Standard, EDI water, RO ^{1st} water
	FEDE-100	Medium pure water system,100L/h, Standard, EDI water, RO ^{1st} water
	FEDE-150	Medium pure water system,1500L/h, Standard, EDI water, RO ^{1st} water
	FEDE-50UT	Medium pure water system,50L/h, Eliminating bacteria and particle, EDI water, RO ^{1st} water
	FEDE-100UT	Medium pure water system,100L/h, Eliminating bacteria and particle, EDI water, RO ^{1st} water
	FEDE-150UT	Medium pure water system,150L/h, Eliminating bacteria and particle, EDI water, RO ^{1st} water

Ordering Information

Host	Item No	Product description			
	FERS-50	Medium double RO water system,50L/h, Standard, RO ^{2nd} water, RO ^{1st} water			
	FERS-100	Medium double RO water system,100L/h, Standard, RO ^{2nd} water, RO ^{1st} water			
	FERS-150	Medium double RO water system,150L/h, Standard, RO ^{2nd} water, RO ^{1st} water			
	FERS-50UT	Medium double RO water system,50L/h, Eliminating bacteria and particle, RO ^{2nd} water, RO ^{1st} water			
	FERS-100UT	Medium double RO water system,100L/h, Eliminating bacteria and particle, RO ^{2nd} water, RO ^{1st} water			
FERS-150UT	Medium double RO water system,150L/h, Eliminating bacteria and particle, RO ^{2nd} water, RO ^{1st} water				
Cartridge	Item No	Product description		Item No	Product description
	HPC104	Pretreatment cartridge D		HPC603	UP cartridge, standard
	HPC366	RO ^{1st} module S660		HPC604	UP cartridge, Low TOC
	HPC393	RO ^{1st} module S912		HPC605	UP cartridge, standard
	HPC399	RO ^{1st} module S918		HPC606	UP cartridge, Low TOC
	HPC397	RO ^{1st} module S925		HPC710	Air filter for tank
	HPC395	RO ^{1st} module F905		HPC703	185&254nm double wavelength UV lamp
	HPC391	RO ^{1st} module F910		HPC704	254nm UV lamp
	HPC396	RO ^{1st} module F915		HPC709	UF ultrafiltration module
	HPC465	RO ^{2nd} module D650		HPC801	TF terminal microfilter
	HPC491	RO ^{2nd} module D910		HPC802	TF terminal microfilter
	HPC496	RO ^{2nd} module D915		HPC803	TF terminal microfilter
	HPC502	DI cartridge		HPC804	TF terminal microfilter
	HPC503	DI cartridge		HPC810	UF terminal ultrafilter
Accessory	Item No	Product description		Item No	Product description
	DISP2001	HiDis dispenser arm (independent), equipped with 2M connection kit		PWA7502	External leak sensor
	PWA7201	Automatic water softener (salt required)		PWSL1000	Pure water supplying module (BA)
	PWA7020	Pretreatment filter for source water		PWSL2000	Pure water storage and supplying module (PR)
	PWA7021	PP cartridge for pretreatment filter		PWSM0500	Pure water storage and supplying module (FC-500L)
	PWA7022	RS cartridge for pretreatment filter		PWSM1000	Pure water storage and supplying module (FC-1000L)
	PWA7501	Foot switch		PWSM2000	Pure water storage and supplying module (FC-2000L)
Service	Item No	Product description			
	HPS51001	1-year extended warranty service (except for consumables)			
	HPS51003	3-year extended warranty service (except for consumables)			
	HPS52001	Verification documents in English			
	HPS53001	Basic verification service			
	HPS59001*	1-year, one-price all-inclusive maintenance agreement, including regular consumables replacement, maintenance and calibration			
HPS59003*	3-year, one-price all-inclusive maintenance agreement, including regular consumables replacement, maintenance and calibration				

*On the basis of mutual confirmation of pure water consumption and feed water quality.



For more product details, please login: www.hpurex.com

Tel: 0086-21-3107 5991

All content in this document (including but not limited to product pictures, Hyperpurex® company LOGO etc.) is owned by Hyperpurex and shall not be used or reproduced without permission.
Hyperpurex shall not be liable for any errors in this document or for any indirect losses arising from the provision, display or use of this document.
The information, instructions and technical specifications in this document are subject to change without notice.

Hyperpurex Instrument Technology (Shanghai) Co., LTD Version: 202503

