

# MULTIPARAMETER METERS



Infitek



Infitek



Infitek



Infitek



Infitek

## Infitek


### Infitek Co., Ltd.

TEL: +86-531-88982330  
FAX: +86-531-88983691  
Website: infitek.com  
Email: info@infitek.com  
Service: support@infitek.com  
Address: Room 201, Building A, No.1 Qianwan 1st Road, Qianhai  
Shenzhen-Hong Kong Cooperation Zone, Shenzhen (Settled in  
Shenzhen Qianhai Business Secretary Co., Ltd.)

US Office

### INFITEK INC.

522W RIVERSIDE AVE STE N,SPOKANE,WA 99201  
EMAIL: INFO@INFITEK.COM

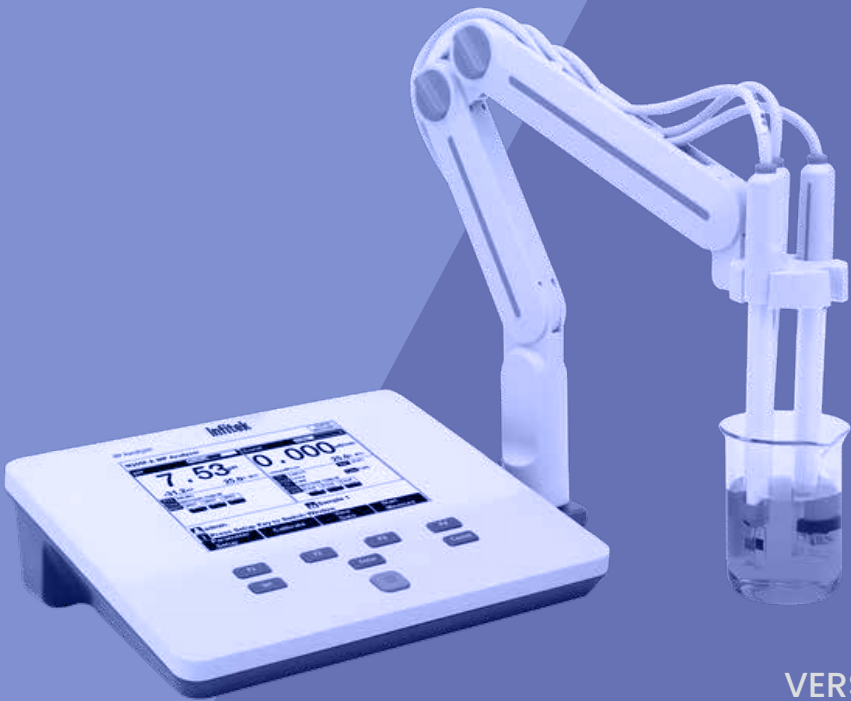

VERSION.2025

# Benchtop Multiparameter Meter

BEP-M300A



## Features

- High resolution LCD display.
- Data Storage 500 sets.
- IP54 waterproof.
- Support for USB and RS-232 communication.



## Description

- The BEP-M300A Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## Features

### General Features

- High resolution LCD display, 5.7 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data Storage 500 sets (GLP-compliant).
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.

### pH

- 1-5 points calibration with Standard Recognition.
- Selectable pH buffer groups, including NIST, DIN, GB.
- Automatic electrode diagnosis with pH slope and offset display.

### Ion

- 1-5 points calibration.
- Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ , PX, etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

### Conductivity

- 1-3 points calibration with Standard Recognition.
- Settable parameters, including cell constant, temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).



### BEP-M300A Multi-parameter Analyzer 1 Set

- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-1VTC Conductivity electrode 1 Piece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standard 1413 $\mu\text{S/cm}$ (50mL) 1 Bottle
- REX-5 electrode stand 1 Set
- European plug model(QD24-0902000-A1)Power adapter DC 9V (Positive inside and negative outside) 1 Piece
- Universal USB cable 1 Piece
- Dust cover (Key function version) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece

Specifications

| Model             |                        | BEP-M300A                                         |
|-------------------|------------------------|---------------------------------------------------|
| Parameters        |                        | pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)  |
| pH                | Range                  | -2.00 to 20.00pH                                  |
|                   | Resolution             | 0.1, 0.01 pH                                      |
|                   | Accuracy               | ±0.01 pH                                          |
|                   | Calibration Points     | Up to 5                                           |
|                   | Standard Customization | Yes                                               |
|                   | Standard Recognition   | NIST, GB and DIN buffers                          |
|                   | Slope Limit            | Yes                                               |
| mV                | Range                  | -2000.0 to 2000.0 mV                              |
|                   | Resolution             | 0.1 mV                                            |
|                   | Accuracy               | ±0.3 mV or ±0.1%                                  |
| pX                | Range                  | -2.00 to 20.00                                    |
|                   | Resolution             | 0.1, 0.01 pX                                      |
|                   | Accuracy               | ±0.01 pX                                          |
|                   | Calibration Points     | Up to 5                                           |
| ISE               | Range                  | 1E-9 to 9.999E9                                   |
|                   | Unit                   | mol/L, mmol/L, g/L, mg/L, µg/L, ppm               |
|                   | Resolution             | Up to 4 significant digits                        |
|                   | Accuracy               | ±0.5%                                             |
|                   | Calibration Points     | Up to 5                                           |
| Conduc-<br>tivity | Range                  | 0.000 µS/cm to 1000 mS/cm                         |
|                   | Resolution             | 0.001 µS/cm minimum, various with range selection |
|                   | Accuracy               | ±1.0% FS                                          |
|                   | Reference Temperature  | 20, 25 °C                                         |
|                   | Calibration Points     | Up to 3                                           |
|                   | Standard Recognition   | 84µS/cm, 1413µS/cm, 12.88mS/cm                    |
| Resistivity       | Range                  | 5.00 Ω·cm~20.00 MΩ·cm                             |
|                   | Resolution             | 0.01 Ω·cm minimum                                 |
|                   | Accuracy               | ±1.0% FS                                          |
| TDS               | Range                  | 0.00ppm~300ppt                                    |
|                   | Resolution             | 0.01mg/L minimum, various with range selection    |
|                   | Accuracy               | ±1.0% FS                                          |
| Salinity          | Range                  | 0~80ppt                                           |
|                   | Resolution             | 0.1ppt                                            |
|                   | Accuracy               | ±2ppt                                             |

| Model                |                               | BEP-M300A                                         |
|----------------------|-------------------------------|---------------------------------------------------|
| Temperature          | Range                         | -5 to 110 °C, 23 to 230 °F                        |
|                      | Unit                          | °C, °F                                            |
|                      | Resolution                    | 0.1                                               |
|                      | Accuracy                      | ±0.2                                              |
| Measurement          | Reading Mode                  | Auto Read (Fast, Medium, Slow), Timed, Continuous |
|                      | Reading Prompts               | Reading, Stable, Locked                           |
|                      | Temp. Compensation            | ATC, MTC                                          |
| Data manage-<br>ment | Data Storage                  | 500 results each                                  |
|                      | GLP Features                  | Yes                                               |
| Inputs               | pH Electrode                  | BNC(Q9)                                           |
|                      | DO with Temp. Probe           | 4-pin aviation connector                          |
|                      | Conductivity with Temp. Probe | 5-pin aviation connector                          |
| Outputs              | USB                           | USB 2.0 flash memory device, PC                   |
|                      | RS 232                        | printer                                           |
| Display options      | Backlight                     | Yes                                               |
|                      | Auto Shutdown                 | 1~60 min, off                                     |
|                      | IP Rating                     | IP54                                              |
|                      | Date and Time                 | Yes                                               |
| General              | Electricity                   | AC Adapter, 100~240V AC input, DC9V output        |
|                      | Dimensions                    | 242×195× 68 mm                                    |
|                      | Weight                        | 900g                                              |



# Benchtop Multiparameter Meter

BEP-M300F



## pH

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including NIST, DIN, GB.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

## Conductivity

- ◆ 1-3 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).


## Ion

- ◆ 1-5 points calibration.
- ◆ Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{PX}$ , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

## DO

- ◆ Support for air-saturated water or zero oxygen calibration.
- ◆ Selectable pressure unit, including kPa, mbar, Torr, Atm.
- ◆ Auto barometric pressure compensation
- ◆ Manual Salinity Factor Correction

## Features

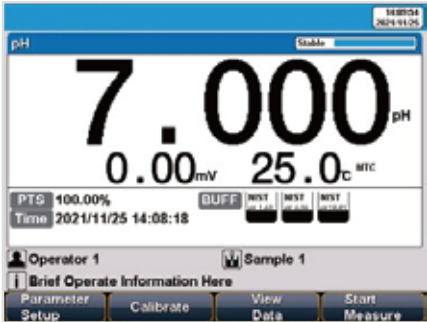

High resolution LCD display, 5.7 inches.  
Smart UI design Best user experience, highly efficient measurement



Support for USB and RS-232C communication.  
Multiple external devices available PC, printer



Automatic identification standard solution  
Fast calibration and more accurate measurement

### Description

The BEP-M300F Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## General Features

IP54 waterproof.

Multi-reading feature allows auto-read, timed-read and continuous-read.

Automatic/Manual temperature compensation ensures accurate results.

Data Storage500 sets (GLP-compliant).

Reset feature automatically resumes all settings back to factory default options.

Auto-hold feature senses and locks the measurement endpoint.



Specifications

| Model        |                        | BEP-M300F                                                         |
|--------------|------------------------|-------------------------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH           | Range                  | -2.00 to 20.00pH                                                  |
|              | Resolution             | 0.1, 0.01 pH                                                      |
|              | Accuracy               | ±0.01 pH                                                          |
|              | Calibration Points     | Up to 5                                                           |
|              | Standard Customization | Yes                                                               |
|              | Standard Recognition   | NIST, GB and DIN buffers                                          |
|              | Slope Limit            | Yes                                                               |
| mV           | Range                  | -2000.0 to 2000.0 mV                                              |
|              | Resolution             | 0.1 mV                                                            |
|              | Accuracy               | ±0.3 mV or ±0.1%                                                  |
| pX           | Range                  | -2.00 to 20.00                                                    |
|              | Resolution             | 0.1, 0.01 pX                                                      |
|              | Accuracy               | ±0.01 pX                                                          |
|              | Calibration Points     | Up to 5                                                           |
| ISE          | Range                  | 1E-9 to 9.999E9                                                   |
|              | Unit                   | mol/L, mmol/L, g/L, mg/L, µg/L, ppm                               |
|              | Resolution             | Up to 4 significant digits                                        |
|              | Accuracy               | ±0.5%                                                             |
|              | Calibration Points     | Up to 5                                                           |
| Conductivity | Range                  | 0.000 µS/cm to 1000 mS/cm                                         |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection                  |
|              | Accuracy               | ±1.0% FS                                                          |
|              | Reference Temperature  | 20, 25 °C                                                         |
|              | Calibration Points     | Up to 3                                                           |
|              | Standard Recognition   | 84µS/cm, 1413µS/cm, 12.88mS/cm                                    |
| Resistivity  | Range                  | 5.00 Ω·cm~20.00 MΩ·cm                                             |
|              | Resolution             | 0.01 Ω·cm minimum                                                 |
|              | Accuracy               | ±1.0% FS                                                          |
| TDS          | Range                  | 0.00ppm~300ppt                                                    |
|              | Resolution             | 0.01mg/L minimum, various with range selection                    |
|              | Accuracy               | ±1.0% FS                                                          |
| Salinity     | Range                  | 0.0~80.0ppt                                                       |
|              | Resolution             | 0.1ppt                                                            |
|              | Accuracy               | ±2ppt                                                             |

| Model                          |                                   | BEP-M300F                                         |
|--------------------------------|-----------------------------------|---------------------------------------------------|
| Dissolved Oxygen Concentration | Range                             | 0.00 to 20.00 ppm                                 |
|                                | Resolution                        | 0.01ppm                                           |
|                                | Accuracy                          | ±0.10ppm                                          |
|                                | Calibration Points                | Air-saturated water or zero point                 |
|                                | Barometric Compensation           | Yes                                               |
| Saturation                     | Manual Salinity Factor Correction | Yes                                               |
|                                | Range                             | (0.0 to 200.0)%                                   |
|                                | Resolution                        | 0.1%                                              |
|                                | Accuracy                          | ±2.0%                                             |
| Temperature                    | Range                             | -5 to 110 °C, 23 to 230 °F                        |
|                                | Unit                              | °C, °F                                            |
|                                | Resolution                        | 0.1                                               |
|                                | Accuracy                          | ±0.2                                              |
| Measurement                    | Reading Mode                      | Auto Read (Fast, Medium, Slow), Timed, Continuous |
|                                | Reading Prompts                   | Reading, Stable, Locked                           |
|                                | Temp. Compensation                | ATC, MTC                                          |
| Data manage-ment               | Data Storage                      | 500 results each                                  |
|                                | GLP Features                      | Yes                                               |
| Inputs                         | pH Electrode                      | BNC(Q9)                                           |
|                                | DO with Temp. Probe               | 4-pin aviation connector                          |
|                                | Conductivity with Temp. Probe     | 5-pin aviation connector                          |
| Outputs                        | USB                               | USB 2.0 flash memory device, PC                   |
|                                | RS 232                            | printer                                           |
| Display options                | Backlight                         | Yes                                               |
|                                | Auto Shutdown                     | 1~60 min, off                                     |
|                                | IP Rating                         | IP54                                              |
|                                | Date and Time                     | Yes                                               |
| General                        | Electricity                       | AC Adapter, 100~240V AC input, DC9V output        |
|                                | Dimensions                        | 242×195× 68 mm                                    |
|                                | Weight                            | 900g                                              |



- BEP-M300F Multi-parameter Analyzer 1 Set**
- E-301-QC 3 in 1 pH composite electrode 1 Piece
  - DJS-IVTC Conductivity electrode 1 Piece
  - DO-958-Q Polarographic DO electrode 1 Piece
  - Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
  - Conductivity standard 1413µS/cm(50mL) 1 Bottle
  - REX-5 electrode stand 1 Set
  - Power adapter DC 9V(European plug model, QD24-0902000-A1) 1 Piece
  - Universal USB cable 1 Piece
  - Dust cover (Key function version) 1 Piece
  - Operation quick guide 1 Piece
  - Instruction manual 1 Piece

# Benchtop Multiparameter Meter

BEP-M500A



### pH

- ◆ 1-6 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including DIN, NIST, USA, MERK, JIS, GB.
- ◆ Support for pH standard group customization.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

### Conductivity

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

### Ion

- ◆ 1-6 points calibration
- ◆ Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

## Features



Color High contrast LCD touch-screen, 7.0 inches.



Support for USB and RS-232 communication. Multiple external devices available  
PC, printer



Automatic identification standard solution  
Fast calibration and more accurate measurement

### Description

The BEP-M500A Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## General Features

IP54 waterproof.

Multi-reading feature allows auto-read, timed-read and continuous-read.

Data Storage 1000 sets for each parameter (GLP-compliant).

Intelligent operation system provides features including user management, calibration management, method management, data management and etc.

Data analysis feature helps user review, compare and recalculate results.

Automatic/Manual temperature compensation ensures accurate results.

Reset feature automatically resumes all settings back to factory default options.

Specifications

| Model        |                        | BEP-M500A                                         |
|--------------|------------------------|---------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)  |
| pH           | Range                  | -2.000 to 20.000 pH                               |
|              | Resolution             | 0.1, 0.01, 0.001pH                                |
|              | Accuracy               | ±0.002 pH                                         |
|              | Calibration Points     | Up to 6                                           |
|              | Standard Customization | Yes                                               |
|              | Calibration Reminder   | Yes                                               |
|              | Standard Recognition   | NIST, GB, DIN, USA, MERK, JIS buffers             |
|              | Slope Limit            | Yes                                               |
| mV/ORP       | Range                  | -2000.00 to 2000.00 mV                            |
|              | Resolution             | 0.1, 0.01 mV                                      |
|              | Accuracy               | ±0.1 mV or ±0.03% of reading whichever is greater |
|              | EH ORP Mode            | Yes                                               |
|              | ORP Calibration Points | 1 custom point (relative mV)                      |
| pX           | Range                  | -2.000 to 20.000 pX                               |
|              | Resolution             | 0.1, 0.01, 0.001 pX                               |
|              | Accuracy               | ±0.002 pX                                         |
|              | Calibration Points     | Up to 6                                           |
| ISE          | Range                  | 1E-9 to 9.999E9                                   |
|              | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L, ppm               |
|              | Resolution             | Up to 4 significant digits                        |
|              | Accuracy               | ±0.3%                                             |
|              | Calibration Points     | Up to 6                                           |
| Conductivity | Range                  | 0.000 µS/cm to 2000 mS/cm                         |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection  |
|              | Accuracy               | ±0.5% FS                                          |
|              | Reference Temperature  | 5, 10, 15, 18, 20, 25 °C                          |
|              | Calibration Points     | Up to 5                                           |
|              | Calibration Reminder   | Yes                                               |
|              | Standard Recognition   | 10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm |
| Resistivity  | Range                  | 5.00 Ω·cm~100.00 MΩ·cm                            |
|              | Resolution             | 0.01 Ω·cm minimum                                 |
|              | Accuracy               | ±0.5% FS                                          |
| TDS          | Range                  | 0.000~1000ppt                                     |
|              | Resolution             | 0.001ppm minimum, changed with range              |
|              | Accuracy               | ±0.5%                                             |

| Model            |                               | BEP-M500A                                         |
|------------------|-------------------------------|---------------------------------------------------|
| Salinity         | Type                          | NaCl                                              |
|                  | Range                         | 0.0~80.0ppt                                       |
|                  | Resolution                    | 0.1ppt                                            |
|                  | Accuracy                      | ±1ppt                                             |
| Temperature      | Range                         | -10 to 135 °C, 14 to 275 °F                       |
|                  | Unit                          | °C, °F                                            |
|                  | Resolution                    | 0.1                                               |
|                  | Accuracy                      | ±0.1                                              |
| Measurement      | Reading Mode                  | Auto Read (Fast, Medium, Slow), Timed, Continuous |
|                  | Reading Prompts               | Reading, Stable, Locked                           |
|                  | Temp. Compensation            | ATC, MTC                                          |
| Data manage-ment | Data Storage                  | 1000 Groups                                       |
|                  | GLP Features                  | Yes                                               |
|                  | Log Management                | Yes                                               |
| Inputs           | pH Electrode                  | BNC(Q9)                                           |
|                  | DO with Temp. Probe           | 4-pin aviation connector                          |
|                  | Conductivity with Temp. Probe | 5-pin aviation connector                          |
| Outputs          | USB                           | USB 2.0 flash memory device, PC                   |
|                  | RS 232                        | printer, scanner, auto-sampler                    |
| Display options  | Backlight                     | Yes                                               |
|                  | Auto Shutdown                 | 1~60 min, off                                     |
|                  | IP Rating                     | IP54                                              |
|                  | Date and Time                 | Yes                                               |
| General          | Electricity                   | AC Adapter, 100~240V AC input, DC24V output       |
|                  | Dimensions                    | 220×195×68 mm                                     |
|                  | Weight                        | 950g                                              |



- BEP-M500A Multi-parameter Analyzer 1 Set**
- E-301-QC 3 in 1 pH composite electrode 1 Piece
  - DJS-IVTC Conductivity electrode 1 Piece
  - Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle
  - Conductivity standard 1413µS/cm(50mL) 1 Bottle
  - REX-6 electrode stand 1 Set
  - GC65-240270-F Power adapter DC24V 1 Piece
  - Dust cover (Touch screen version) 1 Piece
  - European standard power cord(3VTJ2) 1 Piece
  - Universal USB cable 1 Piece
  - Touch pen (Capacitive screen) 1 Piece
  - Operation quick guide 1 Piece
  - Instruction manual 1 Piece




# Benchtop Multiparameter Meter

BEP-M500T



## Features

-  Color High contrast LCD touchscreen.
-  Data Storage 1000 sets for each parameter.
-  IP54 waterproof.
-  Support for USB and RS-232 communication.

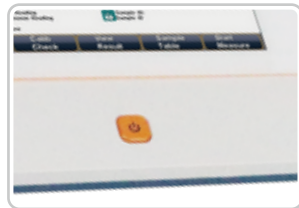

## Description

- The BEP-M500T Benchtop Multi-parameter Analyzer a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## Features

### General Features

- Color High contrast LCD touchscreen, 7.0 inches.
- Intelligent operation system provides features including user management, calibration management, method management, data management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.

### pH

- 1-6 points calibration with Standard Recognition.
- Selectable pH buffer groups, including DIN, NIST, USA, MERK, JIS, GB.
- Support for pH standard group customization.
- Automatic electrode diagnosis with pH slope and offset display.

### Ion

- 1-6 points calibration
- Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

### Conductivity

- 1-5 points calibration with Standard Recognition.
- Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).

### DO

- Support for air-saturated water or zero oxygen calibration.
- Auto barometric pressure compensation.
- Manual Salinity Factor Correction.
- Selectable pressure unit, including  $\text{kPa}$ ,  $\text{mbar}$ ,  $\text{Torr}$ ,  $\text{Atm}$ .



Specifications

| Model         |                        |                                                        | BEP-M500T                                                         |
|---------------|------------------------|--------------------------------------------------------|-------------------------------------------------------------------|
| Parameters    |                        |                                                        | pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH            | Range                  | -2.000 to 20.000 pH                                    |                                                                   |
|               | Resolution             | 0.1, 0.01, 0.001pH                                     |                                                                   |
|               | Accuracy               | ±0.002 pH                                              |                                                                   |
|               | Calibration Points     | Up to 6                                                |                                                                   |
|               | Standard Customization | Yes                                                    |                                                                   |
|               | Calibration Reminder   | Yes                                                    |                                                                   |
|               | Standard Recognition   | NIST, GB, DIN, USA, MERK, JIS buffers                  |                                                                   |
|               | Slope Limit            | Yes                                                    |                                                                   |
| mV/ORP        | Range                  | -2000.00 to 2000.00 mV                                 |                                                                   |
|               | Resolution             | 0.1, 0.01 mV                                           |                                                                   |
|               | Accuracy               | ±0.1 mV or ±0.03% ofreading whichever is greater       |                                                                   |
|               | EH ORP Mode            | Yes                                                    |                                                                   |
|               | ORP Calibration Points | 1 custom point (relative mV)                           |                                                                   |
| pX            | Range                  | -2.000 to 20.000 pX                                    |                                                                   |
|               | Resolution             | 0.1, 0.01, 0.001 pX                                    |                                                                   |
|               | Accuracy               | ±0.002 pX                                              |                                                                   |
|               | Calibration Points     | Up to 6                                                |                                                                   |
| ISE           | Range                  | 1E-9 to 9.999E9                                        |                                                                   |
|               | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L, ppm                    |                                                                   |
|               | Resolution             | Up to 4 significant digits                             |                                                                   |
|               | Accuracy               | ±0.3%                                                  |                                                                   |
|               | Calibration Points     | Up to 6                                                |                                                                   |
| Conduc-tivity | Range                  | 0.000 µS/cm to 2000 mS/cm                              |                                                                   |
|               | Resolution             | 0.001 µS/cm minimum,various with range selection       |                                                                   |
|               | Accuracy               | ±0.5% FS                                               |                                                                   |
|               | Reference Temperature  | 5, 10, 15, 18, 20, 25 °C                               |                                                                   |
|               | Calibration Points     | Up to 5                                                |                                                                   |
|               | Calibration Reminder   | Yes                                                    |                                                                   |
|               | Standard Recognition   | 10 µS/cm, 84 µS/cm, 500 µS/cm, 1413 µS/cm, 12.88 mS/cm |                                                                   |
| Resistivity   | Range                  | 5.00 Ω·cm~100.00 MΩ·cm                                 |                                                                   |
|               | Resolution             | 0.01 Ω·cm minimum                                      |                                                                   |
|               | Accuracy               | ±0.5% FS                                               |                                                                   |
| TDS           | Range                  | 0.000~1000 ppt                                         |                                                                   |
|               | Resolution             | 0.001 ppm minimum, changed with range                  |                                                                   |
|               | Accuracy               | ±0.5%FS                                                |                                                                   |
| Salinity      | Type                   | NaCl                                                   |                                                                   |
|               | Range                  | 0.0~80.0ppt                                            |                                                                   |
|               | Resolution             | 0.1ppt                                                 |                                                                   |
|               | Accuracy               | ±1ppt                                                  |                                                                   |

| Model            |                                   |                                                      | BEP-M500T |
|------------------|-----------------------------------|------------------------------------------------------|-----------|
| Dissolved        | Range                             | 0.00 to 99.99 ppm                                    |           |
| Oxygen           | Resolution                        | 0.01 ppm                                             |           |
| Concentra-tion - | Accuracy                          | ±0.10ppm, 0.00~50.00ppm      ±0.50ppm, 50.0~99.99ppm |           |
|                  | Calibration Reminder              | Yes                                                  |           |
| Polaro-graphic   | Calibration Points                | Air-saturated water or zero point                    |           |
|                  | Barometric Compensation           | Yes                                                  |           |
|                  | Manual Salinity Factor Correction | Yes                                                  |           |
| % Saturation     | Range                             | (0.0 to 600)%                                        |           |
| -Polaro-graphic  | Resolution                        | 0.10%                                                |           |
|                  | Accuracy                          | ±2.0%                                                |           |
| Tempera-ture     | Range                             | -10 to 135 °C, 14 to 275 °F                          |           |
|                  | Unit                              | °C, °F                                               |           |
|                  | Resolution                        | 0.1                                                  |           |
| Measure-ment     | Accuracy                          | ±0.1                                                 |           |
|                  | Reading Mode                      | Auto Read (Fast, Medium, Slow), Timed, Continuous    |           |
|                  | Reading Prompts                   | Reading, Stable, Locked                              |           |
| Data manage-ment | Temp. Compensation                | ATC, MTC                                             |           |
|                  | Data Storage                      | 1000 Groups                                          |           |
|                  | GLP Features                      | Yes                                                  |           |
| Inputs           | Log Management                    | Yes                                                  |           |
|                  | pH Electrode                      | BNC(Q9)                                              |           |
|                  | DO with Temp. Probe               | 4-pin aviation connector                             |           |
| Outputs          | Conductivity with Temp. Probe     | 5-pin aviation connector                             |           |
|                  | USB                               | USB 2.0 flash memory device, PC                      |           |
|                  | RS 232                            | printer, scanner, auto-sampler                       |           |
| Display options  | Backlight                         | Yes                                                  |           |
|                  | Auto Shutdown                     | 1~60 min, off                                        |           |
|                  | IP Rating                         | IP54                                                 |           |
|                  | Date and Time                     | Yes                                                  |           |
| General          | Electricity                       | AC Adapter, 100~240V AC input, DC24V output          |           |
|                  | Dimensions                        | 220×195×68 mm                                        |           |
|                  | Weight                            | 950g                                                 |           |

- BEP-M500T Multi-parameter Analyzer 1 Set**  
E-301-QC 3 in 1 pH composite electrode 1 Piece  
DJS-1VTC Conductivity electrode 1 Piece  
DO-958-Q Polarographic DO electrode 1 Bottle  
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each  
Conductivity standard 1413µS/cm(50mL) 1 Bottle  
REX-6 electrode stand 1 Set  
GC65-240270-F Power adapter DC24V 1 Piece  
Dust cover (Touch screen version) 1 Piece  
European standard power cord(3VTJ2) 1 Piece  
Universal USB cable 1 Piece  
Touch pen (Capacitive screen) 1 Piece  
Operation quick guide 1 Piece  
Instruction manual 1 Piece



# Benchtop Multiparameter Meter

BEP-M600L



## Features

-  Color High contrast LCD touchscreen.
-  Data Storage 1000 sets for each parameter
-  IP54 waterproof.
-  Support for USB and RS-232 communication.



## Description

- The BEP-M600L Benchtop Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## Features

### General Features

- Color High contrast LCD touchscreen, 7.0 inches.
- Intelligent operation system provides features including user management, calibration management, method management, data management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB and RS-232 communication.
- Reset feature automatically resumes all settings back to factory default options.
- IP54 waterproof.
- Support for three modules, namely pH/ISE module, conductivity module and dissolved oxygen module.
- Select from three modules to customize the four channels your way.
- Four channels can be configured in any combination.

### pH

- 1-6 points calibration with Standard Recognition.
- Selectable pH buffer groups, including NIST, DIN, GB.
- Support for pH standard group customization.
- Automatic electrode diagnosis with pH slope and offset display.

### Ion

- 1-6 points calibration
- Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.
- Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

### Conductivity

- 1-5 points calibration with Standard Recognition.
- Settable parameters, including cell constant, reference temperature ( $5/10/15/18/20/25^\circ\text{C}$ ), temperature compensation coefficient and TDS factor.
- Temperature compensation type (none, linear, pure water).


### DO

- Support for air-saturated water or zero oxygen calibration.
- Auto barometric pressure compensation.
- Manual Salinity Factor Correction.
- Selectable pressure unit, including  $\text{kPa}$ ,  $\text{mbar}$ ,  $\text{Torr}$ ,  $\text{Atm}$ .

Specifications

| Model         |                        | BEP-M600L                                                         |
|---------------|------------------------|-------------------------------------------------------------------|
| Parameters    |                        | pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH            | Range                  | -2.000 to 20.000 pH                                               |
|               | Resolution             | 0.1, 0.01, 0.001pH                                                |
|               | Accuracy               | ±0.002 pH                                                         |
|               | Calibration Points     | Up to 6                                                           |
|               | Standard Customization | Yes                                                               |
|               | Calibration Reminder   | Yes                                                               |
|               | Standard Recognition   | NIST, GB, DIN buffers                                             |
|               | Slope Limit            | Yes                                                               |
| mV/ORP        | Range                  | -2000.00 to 2000.00 mV                                            |
|               | Resolution             | 0.1, 0.01 mV                                                      |
|               | Accuracy               | ±0.03% FS                                                         |
|               | EH ORP Mode            | Yes                                                               |
|               | ORP Calibration Points | 1~2 custom point (relative mV)                                    |
| pX            | Range                  | -2.000 to 20.000 pX                                               |
|               | Resolution             | 0.1, 0.01, 0.001 pX                                               |
|               | Accuracy               | ±0.002 pX                                                         |
|               | Calibration Points     | Up to 6                                                           |
| ISE           | Range                  | 1.000E-9 to 9.999E9                                               |
|               | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L                                    |
|               | Resolution             | Up to 4 significant digits                                        |
|               | Accuracy               | ±0.3%                                                             |
|               | Calibration Points     | Up to 6                                                           |
| Conduc-tivity | Range                  | 0.000 µS/cm to 3000 mS/cm                                         |
|               | Resolution             | 0.001 µS/cm minimum, various with range selection                 |
|               | Accuracy               | ±0.5% FS                                                          |
|               | Reference Temperature  | 5, 10, 15, 18, 20, 25 °C                                          |
|               | Calibration Points     | Up to 5                                                           |
|               | Standard Recognition   | 146.5µS/cm, 1408µS/cm, 12.85mS/cm, 111.3mS/cm                     |
| Resistivity   | Range                  | 5.00 Ω·cm~100.00 MΩ·cm                                            |
|               | Resolution             | 0.01 Ω·cm minimum                                                 |
|               | Accuracy               | ±0.5% FS                                                          |
| TDS           | Range                  | 0.000 mg/L~1000 g/L                                               |
|               | Resolution             | 0.001ppm minimum, changed with range                              |
|               | Accuracy               | ±0.5%FS                                                           |
| Salinity      | Type                   | NaCl%                                                             |
|               | Range                  | (0.00~8.00)%                                                      |
|               | Resolution             | 0.01%                                                             |
|               | Accuracy               | ±0.1%                                                             |

| Model            |                                   | BEP-M600L                                         |
|------------------|-----------------------------------|---------------------------------------------------|
| Dissolved        | Range                             | 0.00 to 99.99 mg/L                                |
| Oxygen           | Resolution                        | 0.01 mg/L                                         |
| Concentra-tion - | Accuracy                          | ±0.1 mg/L                                         |
|                  | Calibration Reminder              | Yes                                               |
| Polaro-graphic   | Calibration Points                | Air-saturated water or zero point                 |
|                  | Barometric Compensation           | Yes                                               |
|                  | Manual Salinity Factor Correction | Yes                                               |
| % Saturation     | Range                             | (0.0 to 600)%                                     |
| -Polaro-graphic  | Resolution                        | 0.10%                                             |
|                  | Accuracy                          | ±2.0%                                             |
| Tempera-ture     | Range                             | -5 to 130 °C                                      |
|                  | Unit                              | °C                                                |
|                  | Resolution                        | 0.1                                               |
| Measure-ment     | Accuracy                          | ±0.1                                              |
|                  | Reading Mode                      | Auto Read (Fast, Medium, Slow), Timed, Continuous |
|                  | Reading Prompts                   | Reading, Stable, Locked                           |
| Data manage-ment | Temp. Compensation                | ATC, MTC                                          |
|                  | Data Storage                      | 1000 Groups                                       |
|                  | GLP Features                      | Yes                                               |
| Inputs           | Log Management                    | Yes                                               |
|                  | pH Electrode                      | BNC(Q9)                                           |
|                  | DO with Temp. Probe               | 6-pin Mini-DIN                                    |
| Outputs          | Conductivity with Temp. Probe     | 8-pin Mini-DIN                                    |
|                  | USB                               | USB 2.0 flash memory device, printer, PC          |
| Display options  | RS 232                            | Yes                                               |
|                  | Auto Shutdown                     | 1~60 min, off                                     |
|                  | IP Rating                         | IP54                                              |
|                  | Date and Time                     | Yes                                               |
| General          | Electricity                       | AC Adapter, 100-240V AC input, DC24V output       |
|                  | Dimensions                        | 280×280×130 mm                                    |
|                  | Weight                            | 2500g                                             |

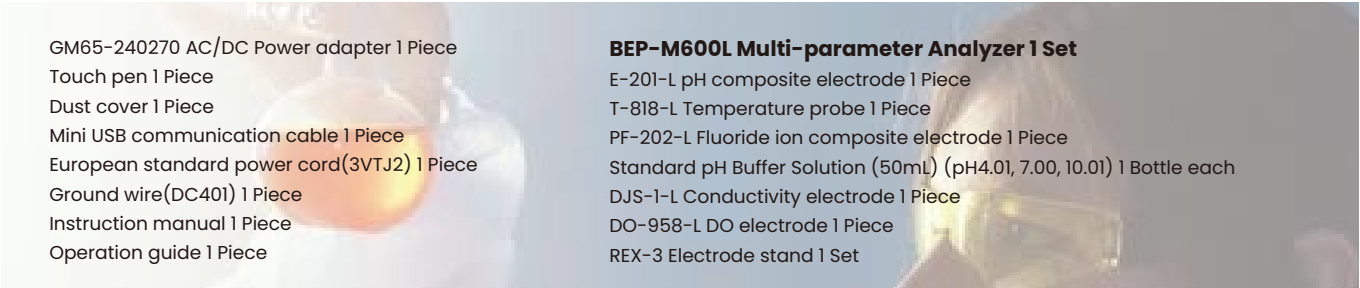

GM65-240270 AC/DC Power adapter 1 Piece

Touch pen 1 Piece

Dust cover 1 Piece

Mini USB communication cable 1 Piece

European standard power cord(3VTJ2) 1 Piece

Ground wire(DC401) 1 Piece

Instruction manual 1 Piece

Operation guide 1 Piece

**BEP-M600L Multi-parameter Analyzer 1 Set**

E-201-L pH composite electrode 1 Piece

T-818-L Temperature probe 1 Piece

PF-202-L Fluoride ion composite electrode 1 Piece

Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each

DJS-1-L Conductivity electrode 1 Piece

DO-958-L DO electrode 1 Piece

REX-3 Electrode stand 1 Set

# Benchtop Multiparameter Meter

BEP-M900



## Measurement Parameters

**BEP-M900:**  
pH, mV, relative mV, ion concentration, conductivity, TDS, salinity, resistivity, DO, temperature.

## General Features



Auto-Read function senses and locks the measurement endpoint



Automatic temperature compensation ensures accurate readings over the entire range



Calibration due alarm prompts the user to calibrate the meter regularly



Setup menu allows setting the number of calibration points, resolution, stability criteria, temperature unit, date and time, etc.



Reset function automatically resumes all settings back to the factory defaults



Expanded memory stores or recalls up to 500 data sets  
USB communication interface for data transfer and timed interval readings

## Features

### pH:

Multiparameter water quality meter is equipped with a 6.5 inches backlit LCD display  
1 to 5 points calibration with automatic recognition for USA, NIST and DIN buffers  
Automatic electrode diagnosis shows the pH slope and offset

### ORP:

1 point offset calibration allows adjusting the displayed value to a known standard  
Relative and absolute millivolt modes ensure the reliable ORP measurements

### Ion Concentration:

2 to 5 points calibration, including the 8 concentration points can be selected  
Direction concentration readout simplifies the measurement process  
Automatic electrode diagnosis shows the calibration points and slopes  
Selectable concentration units (ppm, mg/L, mol/L, mmol/L) and ionic valency





Conductivity/TDS/Salinity/Resistivity:

1 to 5 points calibration with automatic recognition for conductivity standards

Selectable cell constant, reference temperature, TDS factor, linear and pure water compensations, seawater, and practical salinity measurement modes

Automatic electrode diagnosis shows the calibration points and factors



Dissolved Oxygen:

1 or 2 points calibration using the air-saturated water or zero oxygen solution

Salinity and barometric pressure compensations eliminate the measurement error



Ordering Information

**BEP-M900:**  
pH/conductivity/dissolved oxygen electrodes, temperature probe, pH buffer sachets, conductivity standard solutions, DO electrolyte solution, membrane cap, electrode holder, USB cable, and power adapter.



Specifications

| Model                   |                                                           | BEP-M900 |
|-------------------------|-----------------------------------------------------------|----------|
| pH                      |                                                           |          |
| Range                   | -2.000~20.000pH                                           | ●        |
| Resolution              | 0.001, 0.01, 0.1pH, selectable                            | ●        |
| Accuracy                | ±0.002pH                                                  | ●        |
| Calibration Points      | 1 to 5 points                                             | ●        |
| pH Buffer Options       | USA, NIST, DIN, 2 custom buffers                          | ●        |
| ORP                     |                                                           |          |
| Range                   | ±1999.9mV                                                 | ●        |
| Resolution              | 0.1, 1mV, selectable                                      | ●        |
| Accuracy                | ±0.2mV                                                    | ●        |
| Calibration Points      | 1 point                                                   | ●        |
| Ion                     |                                                           |          |
| Range                   | 0.001~19999 (depending on the range of ISE)               | ●        |
| Resolution              | 0.001, 0.01, 0.1, 1                                       | ●        |
| Accuracy                | ±0.5% F.S. (monovalent), ±1% F.S. (divalent)              | ●        |
| Measurement Units       | ppm, mg/L, mol/L, mmol/L                                  | ●        |
| Calibration Points      | 2 to 5 points (0.001, 0.01, 0.1, 1, 10, 100, 1000, 10000) | ●        |
| Conductivity            |                                                           |          |
| Range                   | 0.01~20.00, 200.0, 2000μS/cm, 20.00, 200.0mS/cm           | ●        |
| Resolution              | 0.001,0.01,0.1,1                                          | ●        |
| Accuracy                | ±0.5% F.S.                                                | ●        |
| Calibration Points      | 1 to 5 points                                             | ●        |
| Calibration Solutions   | 10μS/cm, 84μS/cm, 1413μS/cm, 12.88mS/cm, 111.8mS/cm       | ●        |
| Temperature Coefficient | Linear (0.0~10.0%/°C), pure water                         | ●        |
| Reference Temperature   | 20/25°C                                                   | ●        |
| Cell Constant           | K=0.1, 1, 10 or custom                                    | ●        |

Specifications

| Model                          |                                            | BEP-M900 |
|--------------------------------|--------------------------------------------|----------|
| TDS                            |                                            |          |
| Range                          | 0~10.00, 100.0, 1000ppm, 10.00, 200.0ppt   | ●        |
| Resolution                     | 0.01,0.1,1                                 | ●        |
| Accuracy                       | ±1% F.S.                                   | ●        |
| TDS Factor                     | 0.1~1.0 (default 0.5)                      | ●        |
| Salinity                       |                                            |          |
| Range                          | 0.00~42.00psu, 0.00~80.00ppt               | ●        |
| Resolution                     | 0.01                                       | ●        |
| Accuracy                       | ±1% F.S.                                   | ●        |
| Resistivity                    |                                            |          |
| Range                          | 0.00~20.00MΩ                               | ●        |
| Resolution                     | 0.01                                       | ●        |
| Accuracy                       | ±1% F.S.                                   | ●        |
| DO                             |                                            |          |
| Range                          | 0.00~20.00mg/L, 0.0~200.0% saturation      | ●        |
| Resolution                     | 0.01mg/L, 0.1%                             | ●        |
| Accuracy                       | ±0.2mg/L, ±2.0%                            | ●        |
| Calibration Points             | 1or 2 points                               | ●        |
| Barometric Pressure Correction | 60.0~112.5kPa/450~850mmHg, manual          | ●        |
| Salinity Correction            | 0.0~50.0g/L, manual                        | ●        |
| General Spec.                  |                                            |          |
| Temperature Compensation       | 0~100°C/32~212°F, manual or automatic      | ●        |
| Memory                         | 500 data sets, USB communication interface | ●        |
| Power Requirements             | DC5V power adapter                         | ●        |
| Dimensions and Weight          | 210(L)×188(W)×60(H)mm,1.5kg                | ●        |

# Portable Multiparameter Meter

BEP-M310A



## pH

- ◆ 1-5 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including NIST, DIN, GB.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

## Conductivity

- ◆ 1-3 points calibration with Standard Recognition.
- ◆ Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

## Ion

- ◆ 1-5 points calibration.
- ◆ Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.
- ◆ Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

## Features



LCD display screen, 3.5 inches.



Support for USB communication. Multiple external devices available  
PC, printer



Tilting stand  
Support better reading and operation



Wristbands Safety hold

## Description

The BEP-M310A Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## General Features



IP65 waterproof.



Data capacity of up to 500 sets for each parameter (GLP-compliant).



Multi-reading feature allows auto-read, timed-read and continuous-read.



Auto-power off feature effectively extends the battery service life.  
Reset feature automatically resumes all settings back to factory default options.



Automatic/Manual temperature compensation ensures accurate results.  
Auto-hold feature senses and locks the measurement endpoint.

Specifications

| Model        |                        | BEP-M310A                                        |
|--------------|------------------------|--------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.) |
| pH           | Range                  | -2.00 to 20.00 pH                                |
|              | Resolution             | 0.1, 0.01 pH                                     |
|              | Accuracy               | ±0.01 pH                                         |
|              | Calibration Points     | Up to 5                                          |
|              | Standard Customization | Yes                                              |
|              | Standard Recognition   | NIST, GB and DIN buffers                         |
|              | Slope Limit            | Yes                                              |
| mV           | Range                  | -2000.0 to 2000.0 mV                             |
|              | Resolution             | 0.1                                              |
|              | Accuracy               | ±0.3 mV or ±0.1%                                 |
| pX           | Range                  | -2.00 to 20.00                                   |
|              | Resolution             | 0.1, 0.01 pX                                     |
|              | Accuracy               | ±0.01 pX                                         |
|              | Calibration Points     | Up to 5                                          |
| ISE          | Range                  | 1E-9 to 9.999E9                                  |
|              | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L                   |
|              | Resolution             | Up to 4 significant digits                       |
|              | Accuracy               | ±0.5%                                            |
|              | Calibration Points     | Up to 5                                          |
| Conductivity | Range                  | 0.000 µS/cm to 500 mS/cm                         |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection |
|              | Accuracy               | ±1.0% FS                                         |
|              | Reference Temperature  | 20, 25 °C                                        |
|              | Calibration Points     | Up to 3                                          |
|              | Standard Recognition   | 84µS/cm, 1413µS/cm, 12.88mS/cm                   |
| Resistivity  | Range                  | 5.00 Ω·cm~20.00 MΩ·cm                            |
|              | Resolution             | 0.01 Ω·cm minimum                                |
|              | Accuracy               | ±1.0% FS                                         |
| TDS          | Range                  | 0.00 ppm~300 ppt                                 |
|              | Resolution             | 0.01ppm minimum, various with range selection    |
|              | Accuracy               | ±1.0%FS                                          |
| Salinity     | Range                  | 0.0~80.0 ppt                                     |
|              | Resolution             | 0.1ppt                                           |
|              | Accuracy               | ±2ppt                                            |

| Model            |                               | BEP-M310A                                                                |
|------------------|-------------------------------|--------------------------------------------------------------------------|
| Temperature      | Range                         | -5 to 110 °C, 23 to 230 °F                                               |
|                  | Unit                          | °C, °F                                                                   |
|                  | Resolution                    | 0.1                                                                      |
|                  | Relative Accuracy             | ±0.2                                                                     |
| Measurement      | Reading Mode                  | Auto Read (Fast, Medium, Slow), Timed, Continuous                        |
|                  | Reading Prompts               | Reading, Stable, Locked                                                  |
|                  | Temp. Compensation            | ATC, MTC                                                                 |
| Data manage-ment | Data Storage                  | 500 results each                                                         |
|                  | GLP Features                  | Yes                                                                      |
| Inputs           | pH Electrode                  | BNC(Q9)                                                                  |
| Inputs           | Conductivity with Temp. Probe | 5-pin aviation connector                                                 |
| Outputs          | USB                           | PC, printer                                                              |
| Display options  | Backlight                     | Yes                                                                      |
|                  | Auto Shutdown                 | 300, 600, 1200, 1800, 3600 sec., off                                     |
|                  | IP Rating                     | IP65                                                                     |
|                  | Date and Time                 | Yes                                                                      |
| General          | Electricity                   | Rechargeable Lithium battery, AC Adapter, 100-240V AC input, DC5V output |
|                  | Dimensions                    | W80×D255×H35 mm                                                          |
|                  | Package Dimension (W*D*H)(mm) | W360×D300×H180                                                           |
|                  | N.W./G.W.(kg)                 | 0.4kg/2.8kg                                                              |

BEP-M310A Multi-parameter Analyzer 1 Set

MicroUSB to DB25 printing cable 1 Piece  
Silicone rubber case (Installed on the host) 1 Piece  
Wristband assembly(Small) 1 Piece  
Carrying case(Large) 1 Piece  
Operation quick guide 1 Piece  
Instruction manual 1 Piece



E-301-QC 3 in 1 pH composite electrode 1 Piece  
DJS-IVTC Conductivity electrode 1P iece  
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each  
Conductivity standards 1413µS/cm(50mL) 1 Bottle  
Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece  
Electrode holder 2 Pieces  
MicroUSB cable 1 Piece





# Portable Multiparameter Meter

BEP-M310F



The BEP-M310F Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## Features

- LCD display screen, 3.5 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data capacity of up to 500 sets for each parameter (GLP-compliant).
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.

## Features

- ◆ **pH**  
1-5 points calibration with Standard Recognition.  
Selectable pH buffer groups, including NIST, DIN, GB.  
Automatic electrode diagnosis with pH slope and offset display.
- ◆ **Ion**  
1-5 points calibration.  
Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.  
Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.  
Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.
- ◆ **Conductivity**  
1-3 points calibration with Standard Recognition.  
Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.  
Temperature compensation type (none, linear, pure water).
- ◆ **DO**  
Support for air-saturated water or zero oxygen calibration.  
Auto barometric pressure compensation  
Manual Salinity Factor Correction  
Selectable pressure unit, including  $\text{kPa}$ ,  $\text{mbar}$ ,  $\text{Torr}$ ,  $\text{Atm}$ .



LCD display screen.  
Data capacity of up to 500 sets for each parameter.  
Support for USB communication.  
IP65 waterproof.

Specifications

| Model        |                        | BEP-M310F                                                         |
|--------------|------------------------|-------------------------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH           | Range                  | -2.00 to 20.00pH                                                  |
|              | Resolution             | 0.1, 0.01 pH                                                      |
|              | Accuracy               | ±0.01 pH                                                          |
|              | Calibration Points     | Up to 5                                                           |
|              | Standard Customization | Yes                                                               |
|              | Standard Recognition   | NIST, GB and DIN buffers                                          |
|              | Slope Limit            | Yes                                                               |
| mV           | Range                  | -2000.0 to 2000.0 mV                                              |
|              | Resolution             | 0.1                                                               |
|              | Accuracy               | ±0.3 mV or ±0.1%                                                  |
| pX           | Range                  | -2.00 to 20.00                                                    |
|              | Resolution             | 0.1, 0.01 pX                                                      |
|              | Accuracy               | ±0.01 pX                                                          |
|              | Calibration Points     | Up to 5                                                           |
| ISE          | Range                  | 1E-9 to 9.999E9                                                   |
|              | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L                                    |
|              | Resolution             | Up to 4 significant digits                                        |
|              | Accuracy               | ±0.5%                                                             |
|              | Calibration Points     | Up to 5                                                           |
| Conductivity | Range                  | 0.000 µS/cm to 500 mS/cm                                          |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection                  |
|              | Accuracy               | ±1.0% FS                                                          |
|              | Reference Temperature  | 20, 25 °C                                                         |
|              | Calibration Points     | Up to 3                                                           |
|              | Standard Recognition   | 84µS/cm, 1413µS/cm, 12.88mS/cm                                    |
| Resistivity  | Range                  | 5.00 Ω·cm~20.00 MΩ·cm                                             |
|              | Resolution             | 0.01 Ω·cm minimum                                                 |
|              | Accuracy               | ±1.0% FS                                                          |
| TDS          | Range                  | 0.00 ppm~300 ppt                                                  |
|              | Resolution             | 0.01ppm minimum, various with range selection                     |
|              | Accuracy               | ±1.0%FS                                                           |

BEP-M310F Multi-parameter Analyzer 1 Set

- Electrode holder 2 Pieces
- MicroUSB cable 1 Piece
- MicroUSB to DB25 printing cable 1 Piece
- Silicone rubber case (Installed on the host) 1 Piece
- Wristband assembly(Small) 1 Piece
- Carrying case(Large) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece

- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-1VTC Conductivity electrode 1 Piece
- DO-958-Q Polarographic DO electrode 1 Piece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standards 1413µS/cm(50mL) 1 Bottle
- Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece

Specifications

| Model                               |                                   | BEP-M310F                                                                |
|-------------------------------------|-----------------------------------|--------------------------------------------------------------------------|
| Salinity                            | Range                             | 0.0~80.0 ppt                                                             |
|                                     | Resolution                        | 0.1ppt                                                                   |
|                                     | Accuracy                          | ±2ppt                                                                    |
| Dissolved Oxygen Con-<br>centration | Range                             | 0.00 to 50.00ppm                                                         |
|                                     | Resolution                        | 0.01ppm                                                                  |
|                                     | Accuracy                          | ±0.1ppm                                                                  |
|                                     | Calibration Points                | Air-saturated water or zero point                                        |
|                                     | Barometric Compensation           | Yes                                                                      |
|                                     | Manual Salinity Factor Correction | Yes                                                                      |
| Saturation                          | Range                             | (0.0 to 300.0)%                                                          |
|                                     | Resolution                        | 0.1%                                                                     |
|                                     | Accuracy                          | ±10.0%                                                                   |
| Temperature                         | Range                             | -5 to 110 °C, 23 to 230 °F                                               |
|                                     | Unit                              | °C, °F                                                                   |
|                                     | Resolution                        | 0.1                                                                      |
|                                     | Relative Accuracy                 | ±0.2                                                                     |
| Measurement                         | Reading Mode                      | Auto Read (Fast, Medium, Slow), Timed, Continuous                        |
|                                     | Reading Prompts                   | Reading, Stable, Locked                                                  |
|                                     | Temp. Compensation                | ATC, MTC                                                                 |
| Data manage-<br>ment                | Data Storage                      | 500 results each                                                         |
|                                     | GLP Features                      | Yes                                                                      |
| Inputs                              | pH Electrode                      | BNC(Q9)                                                                  |
|                                     | DO with Temp. Probe               | 4-pin aviation connector                                                 |
|                                     | Conductivity with Temp. Probe     | 5-pin aviation connector                                                 |
| Outputs                             | USB                               | PC, printer                                                              |
| Display options                     | Backlight                         | Yes                                                                      |
|                                     | Auto Shutdown                     | 300, 600, 1200, 1800, 3600 sec., off                                     |
|                                     | IP Rating                         | IP65                                                                     |
|                                     | Date and Time                     | Yes                                                                      |
| General                             | Electricity                       | Rechargeable Lithium battery, AC Adapter, 100-240V AC input, DC5V output |
|                                     | Dimensions                        | 80×255×35 mm                                                             |
|                                     | Weight                            | 400g                                                                     |



# Portable Multiparameter Meter

## BEP-M310T



The BEP-M310T Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

### Features

- 1.1LCD display screen,4.3 inches.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data capacity of up to 500 sets for each parameter (GLP-compliant).
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.

### Features

- **pH**  
1-5 points calibration with Standard Recognition.  
Selectable pH buffer groups, including NIST, DIN, GB.  
Automatic electrode diagnosis with pH slope and offset display.
- **Ion**  
1-5 points calibration.  
Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.  
Measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.  
Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.
- **Conductivity**  
1-3 points calibration with Standard Recognition.  
Settable parameters, including cell constant, temperature compensation coefficient, and TDS factor.  
Temperature compensation type (none, linear, pure water).
- **DO**  
Support for air-saturated water or zero oxygen calibration.  
Auto barometric pressure compensation  
Manual Salinity Factor Correction  
Selectable pressure unit, including kPa, mbar, Torr, Atm.



LCD display screen.  
Data capacity of up to 500 sets for each parameter.  
Support for USB communication.  
IP65 waterproof.

Specifications

| Model        |                        | BEP-M310T                                                     |
|--------------|------------------------|---------------------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/DO/Temp. (mV/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH           | Range                  | -2.00 to 20.00pH                                              |
|              | Resolution             | 0.1, 0.01 pH                                                  |
|              | Accuracy               | ±0.01 pH                                                      |
|              | Calibration Points     | Up to 5                                                       |
|              | Standard Customization | Yes                                                           |
|              | Calibration Reminder   | Yes                                                           |
|              | Standard Recognition   | NIST, GB and DIN buffers                                      |
|              | Slope Limit            | Yes                                                           |
|              |                        |                                                               |
| mV           | Range                  | -2000.0 to 2000.0 mV                                          |
|              | Resolution             | 0.1 mV                                                        |
|              | Accuracy               | ±0.3 mV or ±0.1%                                              |
| pX           | Range                  | -2.00 to 20.00 pX                                             |
|              | Resolution             | 0.1, 0.01 pX                                                  |
|              | Accuracy               | ±0.01 pX(±0.01 pX)                                            |
|              | Calibration Points     | Up to 5                                                       |
| ISE          | Range                  | 1E-9 to 9.999E9                                               |
|              | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L                                |
|              | Resolution             | Up to 4 significant digits                                    |
|              | Accuracy               | ±0.5%                                                         |
|              | Calibration Points     | Up to 5                                                       |
| Conductivity | Range                  | 0.000 µS/cm to 1000 mS/cm                                     |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection              |
|              | Accuracy               | ±1.0% FS                                                      |
|              | Reference Temperature  | 20, 25 °C                                                     |
|              | Calibration Points     | Up to 3                                                       |
|              | Standard Recognition   | 84µS/cm, 1413µS/cm, 12.88mS/cm                                |
| Resistivity  | Range                  | 5.00 Ω ·cm~20.00 MΩ ·cm                                       |
|              | Resolution             | 0.01Ω ·cm minimum,                                            |
|              | Accuracy               | ±1.0% FS                                                      |
| TDS          | Range                  | 0.00ppm~300ppt                                                |
|              | Resolution             | 0.01mg/L minimum, 3 significant digits                        |
|              | Accuracy               | ±1.0%                                                         |

|                                                                                                                                                                                                                                                                                                                                                           |                                                                                                                                                                                                                                                                                                                                                                        |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>BEP-M310T Multi-parameter Analyzer 1 Set</b><br>MicroUSB cable 1 Piece<br>MicroUSB to DB25 printing cable 1 Piece<br>Touch pen (Capacitive screen) 1 Piece<br>Silicone rubber case (Installed on the host) 1 Piece<br>Wristband assembly(Large) 1 Piece<br>Carrying case(Large) 1 Piece<br>Operation quick guide 1 Piece<br>Instruction manual 1 Piece | E-301-QC 3 in 1 pH composite electrode 1 Piece<br>DJS-IVTC Conductivity electrode 1 Piece<br>DO-958-Q Polarographic DO electrode 1 Piece<br>Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each<br>Conductivity standards 1413µS/cm(50mL) 1 Bottle<br>Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece<br>Electrode holder 2 Pieces |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

| Model                               |                                   | BEP-M310T                                                                   |
|-------------------------------------|-----------------------------------|-----------------------------------------------------------------------------|
| Salinity                            | Range                             | 0.0~80.0ppt                                                                 |
|                                     | Resolution                        | 0.1ppt                                                                      |
|                                     | Accuracy                          | ±2ppt                                                                       |
| Dissolved Oxygen Con-<br>centration | Sensor Type                       | Polarographic                                                               |
|                                     | Range                             | 0.00 to 50.00 ppm                                                           |
|                                     | Resolution                        | 0.01ppm                                                                     |
|                                     | Accuracy                          | ±0.1ppm                                                                     |
|                                     | Calibration Points                | Air-saturated water or zero point                                           |
|                                     | Barometric Compensation           | Yes                                                                         |
|                                     | Manual Salinity Factor Correction | Yes                                                                         |
| Saturation                          | Range                             | (0.0 to 300.0)%                                                             |
|                                     | Resolution                        | 0.10%                                                                       |
|                                     | Accuracy                          | ±2.0%                                                                       |
| Temperature                         | Range                             | -5 to 110 °C, 23.0 to 230 °F                                                |
|                                     | Unit                              | °C, °F                                                                      |
|                                     | Resolution                        | 0.1                                                                         |
|                                     | Relative Accuracy                 | ±0.2                                                                        |
| Measurement                         | Reading Mode                      | Auto-Read (Fast, Medium, Slow), Timed, Continuous                           |
|                                     | Reading Prompts                   | Reading, Stable, Locked                                                     |
|                                     | Temp. Compensation                | ATC, MTC                                                                    |
| Data manage-<br>ment                | Data Storage                      | 1000 Groups                                                                 |
|                                     | GLP Features                      | Yes                                                                         |
| Inputs                              | pH Electrode                      | BNC(Q9)                                                                     |
|                                     | DO with Temp. Probe               | 4-pin aviation connector                                                    |
|                                     | Conductivity with Temp. Probe     | 5-pin aviation connector                                                    |
| Outputs                             | USB                               | PC, printer                                                                 |
|                                     | Backlight                         | Yes                                                                         |
| Display options                     | Auto Shutdown                     | 300, 600, 1200, 1800, 3600sec., off                                         |
|                                     | IP Rating                         | IP65                                                                        |
|                                     | Date and Time                     | Yes                                                                         |
| General                             | Electricity                       | Rechargeable Lithium battery, AC Adapter, 100-240V<br>AC input, DC5V output |
|                                     | Dimensions                        | 90 x 255 x 40 mm                                                            |
|                                     | Weight                            | 500g                                                                        |





# Portable Multiparameter Meter

BEP-M510A



## pH

- ◆ 1-8 points calibration with Standard Recognition.
- ◆ Selectable pH buffer groups, including USA, NIST, DIN, GB, MERK, JIS.
- ◆ User-defined pH buffer or buffer group is supported.
- ◆ Automatic electrode diagnosis with pH slope and offset display.

## Conductivity

- ◆ 1-5 points calibration automatically recognizes the standard solutions.
- ◆ Settable parameters, including cell constant, reference temperature (5/10/15/18/20/25°C), temperature compensation coefficient and TDS factor.
- ◆ Temperature compensation type (none, linear, pure water).

## Ion

- ◆ 1-6 points calibration.
- ◆ Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.
- ◆ Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.
- ◆ Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.

## Features



Color High contrast LCD touchscreen, 4.3 inches.

Support for USB communication. Multiple external devices available  
PC, printer



Support for USB communication. Multiple external devices available  
PC, printer



Tilting stand  
Support better reading and operation



Wristbands Safety hold

## Description

The BEP-M510A Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## General Features



IP65 waterproof.



Intelligent operation system provides features including user management, calibration management, electrode management, method management, data management, Log Management and etc.



Data Storage 1000 sets for each parameter (GLP-compliant).  
Data analysis feature helps user review, compare and recalculate results.



Multi-reading feature allows auto-read, timed-read and continuous-read.



Automatic/Manual temperature compensation ensures accurate results.  
Auto-hold feature senses and locks the measurement endpoint.



Auto-power off feature effectively extends the battery service life.  
Reset feature automatically resumes all settings back to factory default options.

Specifications

| Model        |                        | BEP-M510A                                         |
|--------------|------------------------|---------------------------------------------------|
| Parameters   |                        | pH/EC/ISE/Temp. (mV/ORP/pX/Resistivity/TDS/Sal.)  |
| pH           | Range                  | -2.000 to 20.000 pH                               |
|              | Resolution             | 0.1, 0.01, 0.001pH                                |
|              | Accuracy               | ±0.002 pH                                         |
|              | Calibration Points     | Up to 8                                           |
|              | Standard Customization | Yes                                               |
|              | Calibration Reminder   | Yes                                               |
|              | Standard Recognition   | USA, NIST, MERK, JIS, GB and DIN buffers          |
|              | Slope Limit            | Yes                                               |
| mV/ORP       | Range                  | -2000.00 to 2000.00 mV                            |
|              | Resolution             | 0.1, 0.01 mV                                      |
|              | Accuracy               | ±0.1 mV or ±0.03%                                 |
|              | EH ORP Mode            | Yes                                               |
|              | ORP Calibration Points | 1 custom point (relative mV)                      |
| pX           | Range                  | -2.000 to 20.000 pX                               |
|              | Resolution             | 0.1, 0.01, 0.001 pX                               |
|              | Accuracy               | ±0.002 pX                                         |
|              | Calibration Points     | Up to 6                                           |
| ISE          | Range                  | 1E-9 to 9.999E9                                   |
|              | Units                  | mol/L, mmol/L, g/L, mg/L, µg/L                    |
|              | Resolution             | Up to 4 significant digits                        |
|              | Accuracy               | ±0.3%                                             |
|              | Calibration Points     | Up to 8                                           |
| Conductivity | Range                  | 0.000 µS/cm to 3000 mS/cm                         |
|              | Resolution             | 0.001 µS/cm minimum,various with range selection  |
|              | Accuracy               | ±0.5% FS                                          |
|              | Reference Temperature  | 5, 10, 15, 18, 20, 25 °C                          |
|              | Calibration Points     | Up to 5                                           |
|              | Calibration Reminder   | Yes                                               |
|              | Standard Recognition   | 10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm |
| Resistivity  | Range                  | 5.00 Ω·cm~100.00 MΩ·cm                            |
|              | Resolution             | 0.01 Ω·cm minimum                                 |
|              | Accuracy               | ±0.5% FS                                          |
| TDS          | Range                  | 0.000ppm~1000ppt                                  |
|              | Resolution             | 0.001ppm minimum, changed with range              |
|              | Accuracy               | ±0.5%FS                                           |

| Model            |                               | BEP-M510A                                                                |
|------------------|-------------------------------|--------------------------------------------------------------------------|
| Salinity         | Range                         | 0.0~80.0ppt                                                              |
|                  | Resolution                    | 0.1ppt                                                                   |
|                  | Accuracy                      | ±1ppt                                                                    |
| Temperature      | Range                         | -10 to 135 °C, 14 to 275 °F                                              |
|                  | Unit                          | °C, °F                                                                   |
|                  | Resolution                    | 0.1                                                                      |
|                  | Relative Accuracy             | ±0.1                                                                     |
| Measurement      | Reading Mode                  | Auto Read (Fast, Medium, Slow), Timed, Continuous                        |
|                  | Reading Prompts               | Reading, Stable, Locked                                                  |
|                  | Temp. Compensation            | ATC, MTC                                                                 |
| Data manage-ment | Data Storage                  | 1000 Groups                                                              |
|                  | GLP Features                  | Yes                                                                      |
| Inputs           | pH Electrode                  | BNC(Q9)                                                                  |
|                  | DO with Temp. Probe           | 4-pin aviation connector                                                 |
|                  | Conductivity with Temp. Probe | 5-pin aviation connector                                                 |
| Outputs          | USB                           | USB 2.0 flash memory device, PC, scanner                                 |
| Display options  | Backlight                     | Yes                                                                      |
|                  | Auto Shutdown                 | 300, 600, 1200, 1800, 3600 sec., off                                     |
|                  | IP Rating                     | IP65                                                                     |
|                  | Date and Time                 | Yes                                                                      |
| General          | Electricity                   | Rechargeable Lithium battery, AC Adapter, 100~240V AC input, DC5V output |
|                  | Dimensions                    | 90×255×40 mm                                                             |
|                  | Package Dimension (W*D*H)     | 360×300×180mm                                                            |
|                  | G.W./N.W.                     | 2.6kg/0.5kg                                                              |

BEP-M510A Multi-parameter Analyzer 1 Set

- OTG adapter cable (MicroUSB connector) 1 Piece
- Touch pen (Capacitive screen) 1 Piece
- Silicone rubber case (Installed on the host) 1 Piece
- Wristband assembly(Large) 1 Piece
- Carrying case(Large) 1 Piece
- Operation quick guide 1 Piece
- Instruction manual 1 Piece



- E-301-QC 3 in 1 pH composite electrode 1 Piece
- DJS-IVTC Conductivity electrode 1 Piece
- Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each
- Conductivity standard 1413µS/cm(50mL) 1 Bottle
- Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece
- Electrode holder 2 Pieces
- MicroUSB cable 1 Piece



# Portable Multiparameter Meter

BEP-M510T



The BEP-M510T Portable Multi-parameter Analyzer adopts a multi-purpose measurement method, which is suitable for the fields of health disease control, environmental protection, biological agriculture and forestry, water quality analysis, scientific research and institutions of higher learning.

## Features

- Color High contrast LCD touchscreen, 4.3 inches.
- Intelligent operation system provides features including user management, calibration management, electrode management, method management, data management, Log Management and etc.
- Multi-reading feature allows auto-read, timed-read and continuous-read.
- Automatic/Manual temperature compensation ensures accurate results.
- Auto-hold feature senses and locks the measurement endpoint.
- Data Storage 1000 sets for each parameter (GLP-compliant).
- Data analysis feature helps user review, compare and recalculate results.
- Support for USB communication.
- Auto-power off feature effectively extends the battery service life.
- Reset feature automatically resumes all settings back to factory default options.
- IP65 waterproof.

## Features

- **pH**  
1-8 points calibration with Standard Recognition.  
Selectable pH buffer groups, including USA, NIST, DIN, GB, MERK, JIS.  
User-defined pH buffer or buffer group is supported.  
Automatic electrode diagnosis with pH slope and offset display.
- **Ion**  
1-6 points calibration.  
Selectable measurement unit, including  $\mu\text{g/L}$ ,  $\text{mg/L}$ ,  $\text{g/L}$ ,  $\text{mmol/L}$ ,  $\text{pX}$ ,  $\text{ppm}$ , etc.  
Multi-measurement modes are supported, including Direct Reading mode, Standard Addition mode, Sample Addition mode and GRAN mode.  
Over 10 methods are built-in, including  $\text{F}^-$ ,  $\text{Cl}^-$ ,  $\text{Br}^-$ ,  $\text{I}^-$ ,  $\text{NO}_3^-$ ,  $\text{BF}_4^-$ ,  $\text{NH}_4^+$ ,  $\text{K}^+$ ,  $\text{Na}^+$ ,  $\text{Ca}^{2+}$ ,  $\text{Cu}^{2+}$ ,  $\text{Pb}^{2+}$ ,  $\text{Ag}^+$  and etc., user-defined method is supported.
- **Conductivity**  
1-5 points calibration automatically recognizes the standard solutions.  
Settable parameters, including cell constant, reference temperature ( $5/10/15/18/20/25^\circ\text{C}$ ), temperature compensation coefficient and TDS factor.  
Temperature compensation type (none, linear, pure water).
- **DO**  
Support for air-saturated water or zero oxygen calibration.  
Auto barometric pressure compensation  
Manual Salinity Factor Correction  
Selectable pressure unit, including  $\text{kPa}$ ,  $\text{mbar}$ ,  $\text{Torr}$ ,  $\text{Atm}$ .



Color High contrast LCD touchscreen.  
Data Storage 1000 sets for each parameter.  
Support for USB communication.  
IP65 waterproof.

Specifications

| Model                |                                                   | BEP-M510T                                                         |
|----------------------|---------------------------------------------------|-------------------------------------------------------------------|
| Parameters           |                                                   | pH/EC/ISE/DO/Temp. (mV/ORP/pX/Resistivity/TDS/Sal./DO Saturation) |
| pH                   | Range                                             | -2.000 to 20.000 pH                                               |
|                      | Resolution                                        | 0.1, 0.01, 0.001pH                                                |
|                      | Accuracy                                          | ±0.002 pH                                                         |
|                      | Calibration Points                                | Up to 8                                                           |
|                      | Standard Customization                            | Yes                                                               |
|                      | Calibration Reminder                              | Yes                                                               |
|                      | Standard Recognition                              | USA, NIST, MERK, JIS, GB and DIN buffers                          |
|                      | Slope Limit                                       | Yes                                                               |
| mV/ORP               | Range                                             | -2000.00 to 2000.00 mV                                            |
|                      | Resolution                                        | 0.1, 0.01 mV                                                      |
|                      | Accuracy                                          | ±0.1 mV or ±0.03%                                                 |
|                      | EH ORP Mode                                       | Yes                                                               |
| pX                   | ORP Calibration Points                            | 1 custom point (relative mV)                                      |
|                      | Range                                             | -2.000 to 20.000 pX                                               |
|                      | Resolution                                        | 0.1, 0.01, 0.001 pX                                               |
|                      | Accuracy                                          | ±0.002 pX                                                         |
| Calibration Points   | Up to 6                                           |                                                                   |
| ISE                  | Range                                             | 1E-9 to 9.999E9                                                   |
|                      | Units                                             | mol/L, mmol/L, g/L, mg/L, µg/L                                    |
|                      | Resolution                                        | Up to 4 significant digits                                        |
|                      | Accuracy                                          | ±0.3%                                                             |
| Calibration Points   | Up to 8                                           |                                                                   |
| Conductivity         | Range                                             | 0.000 µS/cm to 3000 mS/cm                                         |
|                      | Resolution                                        | 0.001 µS/cm minimum,various with range selection                  |
|                      | Accuracy                                          | ±0.5% FS                                                          |
|                      | Reference Temperature                             | 5, 10, 15, 18, 20, 25 °C                                          |
|                      | Calibration Points                                | Up to 5                                                           |
|                      | Calibration Reminder                              | Yes                                                               |
| Standard Recognition | 10µS/cm, 84µS/cm, 500µS/cm, 1413µS/cm, 12.88mS/cm |                                                                   |
| Resistivity          | Range                                             | 5.00 Ω·cm~100.00 MΩ·cm                                            |
|                      | Resolution                                        | 0.01 Ω·cm minimum                                                 |
|                      | Accuracy                                          | ±0.5% FS                                                          |
| TDS                  | Range                                             | 0.000ppm~1000ppt                                                  |
|                      | Resolution                                        | 0.001ppm minimum, changed with range                              |
|                      | Accuracy                                          | ±0.5%FS                                                           |

- BEP-M510T Multi-parameter Analyzer 1 Set**  
MicroUSB cable 1 Piece  
OTG adapter cable (MicroUSB connector) 1 Piece  
Touch pen (Capacitive screen) 1 Piece  
Silicone rubber case (Installed on the host) 1 Piece  
Wristband assembly(Large) 1 Piece  
Carrying case(Large) 1 Piece  
Operation quick guide 1 Piece

Instruction manual 1 Piece  
E-301-QC 3 in 1 pH composite electrode 1 Piece  
DJS-IVTC Conductivity electrode 1 Piece  
DO-958-Q Polarographic DO electrode 1 Piece  
Standard pH Buffer Solution (50mL) (pH4.01, 7.00, 10.01) 1 Bottle each  
Conductivity standard 1413µS/cm(50mL) 1 Bottle  
Power adapter DC 5V 2A(European plug model, GAT-0502000) 1 Piece  
Electrode holder 2 Pieces

| Model                               |                                   | BEP-M510T                                                                |
|-------------------------------------|-----------------------------------|--------------------------------------------------------------------------|
| Salinity                            | Range                             | 0.0~80.0ppt                                                              |
|                                     | Resolution                        | 0.1ppt                                                                   |
|                                     | Accuracy                          | ±1ppt                                                                    |
| Dissolved Oxygen Con-<br>centration | Sensor Type                       | Polarographic                                                            |
|                                     | Range                             | 0.00 to 99.99 ppm                                                        |
|                                     | Resolution                        | 0.01ppm                                                                  |
|                                     | Accuracy                          | ±0.1ppm (0.00 to 50.00 ppm);±0.5ppm (50.00 to 99.99 ppm)                 |
|                                     | Calibration Reminder              | Yes                                                                      |
|                                     | Calibration Points                | Air-saturated water or zero point                                        |
|                                     | Barometric Compensation           | Yes                                                                      |
|                                     | Manual Salinity Factor Correction | Yes                                                                      |
| Saturation                          | Range                             | (0.0 to 600.0)%                                                          |
|                                     | Resolution                        | 0.10%                                                                    |
|                                     | Accuracy                          | ±10.0%                                                                   |
| Temperature                         | Range                             | -10 to 135 °C, 14 to 275 °F                                              |
|                                     | Unit                              | °C, °F                                                                   |
|                                     | Resolution                        | 0.1                                                                      |
|                                     | Relative Accuracy                 | ±0.1                                                                     |
| Measurement                         | Reading Mode                      | Auto Read (Fast, Medium, Slow), Timed, Continuous                        |
|                                     | Reading Prompts                   | Reading, Stable, Locked                                                  |
|                                     | Temp. Compensation                | ATC, MTC                                                                 |
| Data manage-<br>ment                | Data Storage                      | 1000 Groups                                                              |
|                                     | GLP Features                      | Yes                                                                      |
| Inputs                              | pH Electrode                      | BNC(Q9)                                                                  |
|                                     | DO with Temp. Probe               | 4-pin aviation connector                                                 |
|                                     | Conductivity with Temp. Probe     | 5-pin aviation connector                                                 |
| Outputs                             | USB                               | USB 2.0 flash memory device, PC, scanner                                 |
| Display options                     | Backlight                         | Yes                                                                      |
|                                     | Auto Shutdown                     | 300, 600, 1200, 1800, 3600 sec., off                                     |
|                                     | IP Rating                         | IP65                                                                     |
|                                     | Date and Time                     | Yes                                                                      |
| General                             | Electricity                       | Rechargeable Lithium battery, AC Adapter, 100~240V AC input, DC5V output |
|                                     | Dimensions                        | 90×255×40 mm                                                             |
|                                     | Weight                            | 500g                                                                     |





# Portable Multiparameter Meter

## BEP-M900P



### Measurement Parameters

**BEP-M900P:** pH, mV, relative mV, ion concentration, conductivity, TDS,salinity, resistivity, DO, temperature.

### General Features

Calibration due alarm prompts user to calibrate the meter regularly

Expanded memory stores or recalls up to 500 data sets

USB communication interface is easy to transfer data to PC

Multi-mode power scheme (batteries, power adapter and USB port) ensures that using

Setup menu allows setting the number of calibration points, resolution, stability criteria, temperature unit, date and time, etc.

Reset function automatically resumes all settings back to the factory defaults

Automatic temperature compensation ensures accurate readings over the entire range  
Auto-Read function senses and locks the measurement endpoint  
Auto-Power Off effectively conserves battery life

### Features

**pH:**  
Multi parameter water quality meter is equipped with a 3.5 inches backlit LCD display  
1 to 5 points calibration with automatic recognition for USA, NIST and DIN buffers  
Automatic electrode diagnosis shows the pH slope and offset

**ORP:**  
1 point offset calibration allows adjusting the displayed value to a known standard  
Relative and absolute millivolt modes ensure the reliable ORP measurements

**Dissolved Oxygen:**  
1 or 2 points calibration using the air-saturated water or zero oxygen solution  
Salinity and barometric pressure compensations eliminate the measurement error

**Ion Concentration:**  
2 to 5 points calibration, including the 8 concentration points can be selected  
Direct ion concentration readout simplifies the measurement process  
Automatic electrode diagnosis shows the calibration points and slopes  
Selectable concentration units (ppm, mg/L, mol/L, mmol/L) and ionic valency

**Conductivity/TDS/Salinity/Resistivity:**  
1 to 5 points calibration with automatic recognition for conductivity standards  
Selectable cell constant, reference temperature, TDS factor, linear and pure water  
compensations, seawater and practical salinity measurement modes  
Automatic electrode diagnosis shows the calibration points and factors



Specifications

| Model                   |                                                           | BEP-M900P |
|-------------------------|-----------------------------------------------------------|-----------|
| pH                      |                                                           |           |
| Range                   | -2.000~20.000pH                                           | ●         |
| Resolution              | 0001,0.01,0.1pH, selectable                               | ●         |
| Accuracy                | 0.002pH                                                   | ●         |
| Calibration Points      | 1 to 5 points                                             | ●         |
| pHBuffer Options        | USA, NIST, DIN, 2 custom buffers                          | ●         |
| ORP                     |                                                           |           |
| Range                   | ±1999.9mV                                                 | ●         |
| Resolution              | 0.1, 1mV, selectable                                      | ●         |
| Accuracy                | ±0.2mV                                                    | ●         |
| Calibration Points      | 1 point                                                   | ●         |
| Ion                     |                                                           |           |
| Range                   | 0.001~19999 (depending on the range of ISE)               | ●         |
| Resolution              | 0.001, 0.01, 0.1, 1                                       | ●         |
| Accuracy                | ±0.5% F.S. (monovalent), ±1% F.S. (divalent)              | ●         |
| Measurement Units       | ppm, mg/L, mol/L, mmol/L                                  | ●         |
| Calibration Paints      | 2 to 5 points (0.001, 0.01, 0.1, 1, 10, 100, 1000, 10000) | ●         |
| Conductivity            |                                                           |           |
| Range                   | 0.01~20.00, 200.0, 2000μS/cm, 20.00, 200.0mS/cm           | ●         |
| Resolution              | 0001,0.01,0.1,1                                           | ●         |
| Accuracy                | ±0.5% F.S.                                                | ●         |
| Calibration Paints      | 1 to 5 points                                             | ●         |
| Calibration Solutions   | 10μS/cm, 84μS/cm, 1413μS/cm, 12.88mS/cm, 111.8mS/cm       | ●         |
| Temperature Coefficient | Linear(0.0~10.0%/°C).pure water                           | ●         |
| Reference Temperature   | 20/25°C                                                   | ●         |
| Cell Constant           | K=0.1, 1, 10 or custom                                    | ●         |

Specifications

| Model                          |                                            | BEP-M900P |
|--------------------------------|--------------------------------------------|-----------|
| TDS                            |                                            |           |
| Range                          | 0~10.00, 100.0, 1000ppm, 10.00, 200.0ppt   | ●         |
| Resolution                     | 001,0.1,1                                  | ●         |
| Accuracy                       | ±1% F.S.                                   | ●         |
| TDS Factor                     | 0.1~1.0 (default 0.5)                      | ●         |
| Salinity                       |                                            |           |
| Range                          | 0.00~42.00psu, 0.00~80.00ppt               | ●         |
| Resolution                     | 0.01                                       | ●         |
| Accuracy                       | ±1% F.S.                                   | ●         |
| Resistivity                    |                                            |           |
| Range                          | 0.00~20.00MΩ                               | ●         |
| Resolution                     | 0.01                                       | ●         |
| Accuracy                       | ±1% F.S.                                   | ●         |
| DO                             |                                            |           |
| Range                          | 0.00~20.00mg/L, 0.0~200.0% saturation      | ●         |
| Resolution                     | 001mg/L,0.1%                               | ●         |
| Accuracy                       | ±0.2mg/L, ±2.0%                            | ●         |
| Calibration Points             | 1 or 2 points                              | ●         |
| Barometric Pressure Correction | 60.0~112.5kPa/450~850mmHg, manual          | ●         |
| Salinity Correction            | 0.0~50.0g/L, manual                        | ●         |
| General Spec.                  |                                            |           |
| Temperature Compensation       | 0~100°C/32~212°F, manual or automatic      | ●         |
| Memory                         | 500 data sets, USB communication interface | ●         |
| Power Requirements             | 3× 1.5V AA batteries or DC5V power adapter | ●         |
| Dimensions and Weight          | W*D*H170× 85 ×30mm, 300g                   | ●         |
| Package Dimension (W*D*H)(mm)  | W*D*H440*320*120mm                         | ●         |
| G.W.(kg)                       | 2.2kg                                      | ●         |